HighBackSling Disposable

SystemRoMedic™

Manual - English



Manual - Svenska

Manual - Norsk

Manual - Danske

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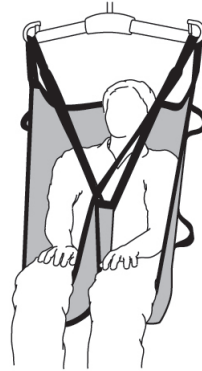
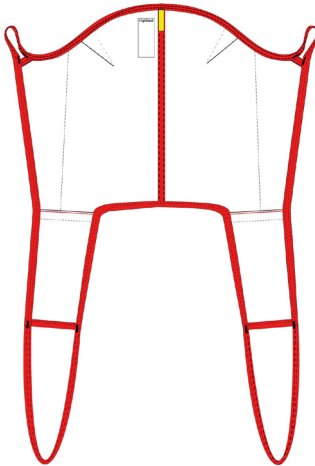
REF

40890004 S
40890006 M
40890007 L
40890008 XL
40890009 XXL

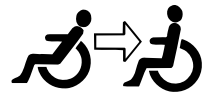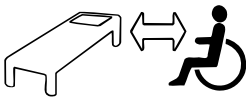SWL:

S-L: 200 kg/440 lbs

XL-XXL: 300 kg/660 lbs

Table of contents

Size	3
English	4-7
Svenska	8-11
Norsk	12-15
Danske	16-19
Deutsch	20-23
Français	24-27

Size

		S/Red	M/Yellow	L/Green	XL/Purple	XXL/White
Rec. User weight	lbs	66-121	110-198	176-243	243-353	330-551
	kg	30-55	50-90	80-110	110-160	150-250
Sling overall height	inch	36	37	38,5	38,5	40
	cm	91	94	98	98	102
Sling seat measurement	inch	14	18	22	24,5	27,5
	cm	35	45	56	62	70
Rec. sling bar measurement	inch	12-21,5	12-23,5	17,5-27,5	17,5-27,5	17,5-27,5
	cm	30-55	30-60	45-70	45-70	45-70

HighBackSling Disposable

SystemRoMedic™

English

HighBackSling Disposable is intended for short-term use. It is easy for the caregiver to apply and it gives the patient a feeling of security. It features leg sections, which means that the patient does not sit on the sling. Instead, the sling provides support from the tailbone, up over the head and under the legs. HighBackSling Disposable is made of a durable, disposable material, which means that it cannot be laundered. Instead, the sling is discarded when it becomes soiled or when the patient no longer needs it.

SystemRoMedic includes a wide range of functional, comfortable, high-quality slings that can be adapted for different types of lifting and for patients with different needs. The slings are available in several different materials and in sizes from XXS to XXL. All models are safe and very easy to use and are rated for lifting patients weighing up to 300 kg. The choice of model and material depends on the patient's needs and the transfer situation.

Functional inspection

Visual and mechanical inspection

Check the condition and function of the sling regularly. The sling must not be used if it has been laundered. Check to ensure that seams and material are free from damage. Check to ensure that fabric is not worn or faded. Apply load to the device and check to ensure that clasps, handles, etc. withstand heavy load. If there are signs of wear, the product must be discarded.

Always read the manual

Always read the manuals for all assistive devices used during a transfer.

Keep the manual where it is accessible to users of the product.

Do not leave the patient unattended during a transfer situation.


To prevent discomfort and the risk of the patient sliding out of the sling, trial fit the sling carefully. First, apply the sling's back loops to the slingbar, and then apply the leg support loops. Use a slingbar that has been tested.

The lowest allowable safe working load always determines the safe working load of the assembled system. Always check the safe working loads for the lift and accessories before use. Contact your dealer if you have any questions.

Application of the sling when the patient is sitting in a wheelchair/chair or in bed


Hold the sling in one hand and place your other hand in the pocket at the lower edge of the sling. This is easier if the patient leans forward or can be helped to lean forward slightly. Guide the sling down behind the patient's back, by sliding it along the back of the chair/wheelchair or along the bed, ideally, to the tailbone. The mid-section of the sling should follow the patient's spine.

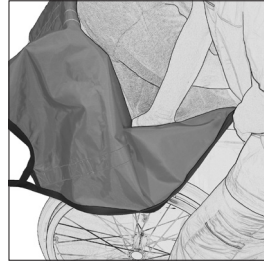Take care to guide the sling down to the patient's hips. Do this by pushing the sling into the space in the back corners of the wheelchair seat.

Then, work the leg sections out along the outside of the patient's thighs. Face the patient and pull the lower edge of the leg support, so that the sling slides into place around the hips. Light pressure against the patient's knees prevents him/her from sliding forward in the chair.

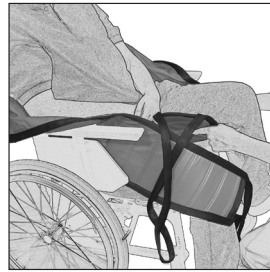After checking to ensure that the leg supports are of equal length on both sides, place them under the patient's thighs. Cross the leg supports and hook them to the slingbar.

Start raising the lift, when the straps become taut, stop and check that everything is ok before proceeding with the lift.

Removal of the sling when the patient is sitting


Fold the leg support and allow it to slide under itself when you pull it out at the outside of the thigh.

Then, pull the sling up behind the patient's back.

Application of the sling when the patient is lying in bed, on a stretcher or on the floor

Turn the patient towards you. Place the sling so that the product label is facing the underlying surface and the sling's mid-section follows the spine. The upper edge of the sling should be high enough so that the sling supports the head and the lower edge is level with the tailbone. Place the sling as far in under the patient as possible. Push the back-section loop under, so that you can reach it once you have turned the patient back onto the sling.

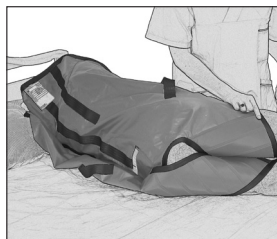To reduce the number of application steps, you can fold the back section up against the patient's back and place the leg supports between the patient's knees when he/she is turned towards you, and then pull them forward when the patient is lying on the sling.

Turn the patient back onto the sling and pull out the head support, back section and leg support.

Cross the leg supports and hook them to the slingbar. Start raising the lift, when the straps become taut, stop and check that everything is ok before proceeding with the lift.


Removal of the sling when the patient is recumbent

Lay the leg supports between the patient's knees.

Turn the patient towards you and fold the sling in under the patient so that the fabric slides towards you when you pull it out after you have turned back the patient.

Alternative application of the sling's leg supports

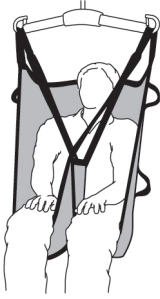In the standard application the leg supports are crossed before they are hooked to the slingbar.

If necessary, the leg supports can also be placed overlapped under both legs, and then passed through the guide loops at the edge of the sling.

For more access to the perinea area of the patient. You can even choose not to cross the leg supports. Note, however, that this entails a risk that the patient might slide out of the sling.

Size information

Slings can often be used for both 2 and 4-point suspension. It is essential to ensure that the width of the slingbar is appropriate for the sling size and the patient's status and capacity.

Use patient-specific slings to ensure safety and hygiene.

Care of the product:

Read the product label.

Do not launder. Discard the product if it is soiled or if there is a need to use another product.

Material:

Non-Woven, polypropylene

CE Medical Device Class I

The product meets the requirements of the Medical Devices Directive 93/42/EEC.

Contact your local distributor if you have any questions about the product and its use.

See www.directhealthcaregroup.com for a complete list of distributors. Always make sure that you have the right version of the manual. The most recent editions of manuals are available for downloading from our website, www.directhealthcaregroup.com.

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Simple solutions for great results

SystemRoMedic™ is the name of Direct Healthcare Group's unique easy transfer concept, the market's widest and most comprehensive range of clever, easy-to-use and safe transfer and lifting aids designed to make life easier, both for the user and for the caregiver. SystemRoMedic™ is a complete solution that provides for the majority of patient transfer or manual handling requirements. From the simplest to the most complex scenarios, from the lightest to the heaviest. The concept encompasses assistive devices for four different categories of transfers:

- Transfer, assistive devices for manual transfers of users between two locations.
- Positioning, assistive devices for manual repositioning of users within the same location.
- Support, assistive devices for mobility support e.g., during sit-to-stand or gait training.
- Lifting, assistive devices for manual and mechanical lifting of users.

Improved work environment, improved quality of care and cost savings

The philosophy behind SystemRoMedic™ is focused on the prevention and reduction of occupational injuries while allowing users to experience a greater sense of independence and dignity. Through a unique combination of training and a complete range of efficient transfer aids, SystemRoMedic™ offers improvement of both work environment and quality of care and, at the same time, achieves significant cost savings.

Always make sure that you have the correct version of the manual

The most recent version of all manuals are available for downloading at/from our website;
www.directhealthcaregroup.com

For questions about the product and its use

Please contact your local Direct Healthcare Group and SystemRoMedic™ representative. A complete list of all our partners with their contact details can be found on our website; www.directhealthcaregroup.com.

Direct Healthcare Group

Advancing Movement & Health

Direct Healthcare Group Sverige AB
Torshamnsgatan 35,
SE-164 40 Kista, Sweden

Tel: +46 (0)8-557 62 200
info@directhealthcaregroup.com
www.directhealthcaregroup.com