# Care instructions

for Direct Healthcare Group's SystemRoMedic® products.

English  
Svenska  
Norsk  
Dansk  
Suomi  
Deutsch  
Nederlands  
Français  
Italiano  
Español

## SystemRoMedic®

SystemRoMedic® is a complete solution for easy transfers. The concept encompasses four product areas that meet most needs in different transfer situations.

**Transfer**, between two locations.

**Positioning**, within the same location.

**Support**, for mobility.


**Lifting**, both manual and mechanical.

**Direct Healthcare  
Group**

Advancing Movement & Health®

## Care instructions

# Using Direct Healthcare Group's SystemRoMedic® products should feel good from start to finish.
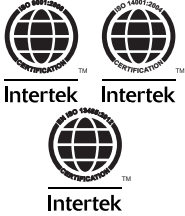
### Quality and environment certification

Direct Healthcare Group is quality and environment certified in accordance with ISO 14001 and ISO 13485. Direct Healthcare Group is continuously and actively working to secure and develop the quality and the environmental thinking within the company.

Direct Healthcare Group strives to ensure that all textile products are Ökotex certified, which ensures that they do not contain any health-hazardous substances in dangerous concentrations.

### Tested and approved products

The assistive devices manufactured by Direct Healthcare Group are of top-quality and fulfill established standards for these types of medical- and healthcare aids. All of our assistive devices are CE-marked and many are patented. All of our assistive devices comply with the requirements of the Medical Devices Regulation for Class 1 products (MDR 2017/745).

### Functional materials

ULF™ (Ultra Low Friction) is the registered trademark for Direct Healthcare Group's sliding material which is used to achieve extremely low friction in, for example, sliding mats.

### Disposable and wipeable materials facilitate infection control efforts

Direct Healthcare Group's SystemRoMedic® range offers an increasing selection of personal transfer aids for short-term use. Instead of laundering, these products are discarded after use, when they get soiled or when the user no longer needs them. This prevents the spread of multi-resistant and other types of contagious bacteria.

Our wipeable transfer aids are made from a strong, safe and comfortable material with a smooth and soft surface and many practical features. The surface of the entire product can be cleaned using disinfectants. In fact, we call this material "wipeable" but these products can actually be dipped in water or used in the shower without getting wet.

### Protective covers

Most sizes of EasySlide can be fitted with a replaceable protective cover, which is made of polyester and can easily be cleaned. Most sizes of EasySlide and RollerSlide can also be fitted with disposable covers, made of bio-degradable polyethylene, which provide for even further improved hygiene and less laundry.

### Laundering and cleaning instructions

Laundry, cleaning and reconditioning of assistive devices is a prerequisite for hygiene and safety as well as for of infection control. General information on how to care for our products can be found hereunder.

For more detailed information, please read the instructions on the product labels. Our products are continuously being developed and improved, which is why we reserve the right to make changes in products without prior notice.

### Textile products

All of our textile manual devices can be washed in 60-80°C/140-176°F water and tumbler-dried at low temperature. Low temperature for tumbler-drying is 45-50°C/104-120°F.

Do not use rinse/fabric softener and avoid tumbler-drying for optimal service life. If disinfection is needed, 70 % ethanol, 45 % isopropanol or similar should be used.

All of our textile lifting slings can be washed in 60-85°C/140-185°F water and tumbler-dried at low temperature.

WendyLett, WendyLean, ReadyLett TopSheet and IncoSheet can be washed in 60-90°C/140-194°F water and tumbler-dried.

### Plastic products

The surface of the product can be cleaned with warm water or a soap solution. If disinfection is needed, 70 % ethanol, 45 % isopropanol or similar should be used.

### Mobile and stationary lifts

If necessary, clean the lift with warm water or a soap solution and check that the castors are free from dirt and hair. Do not use cleaning agents containing phenol or chlorine, as this could damage the aluminium and the plastic materials.

If disinfection is needed, 70 % ethanol, 45 % isopropanol or similar should be used.

## SystemRoMedic®

SystemRoMedic® is a complete solution for easy transfers. The concept encompasses four product areas that meet most needs in different transfer situations.

**Transfer**, between two locations. **Positioning**, within the same location.

**Support**, for mobility. **Lifting**, both manual and mechanical.

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