ShellHighBackSling

SystemRoMedic™

Instructions for use – English

Bruksanvisning – Svenska

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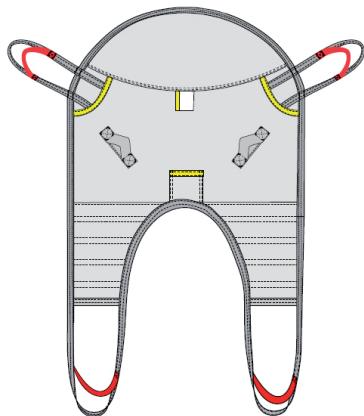
REF

40700002

40700003

40700004

40700006

40700007


40700008

40700009


SWL: 300kg/660lbs

Table of contents

English	3
Svenska	8
Norsk	13
Dansk	18
Deutsch	23
Nederlands	28
Français	33
Italiano	38

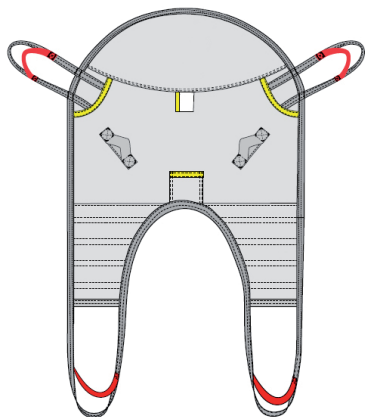SWL: 300kg/660lbs

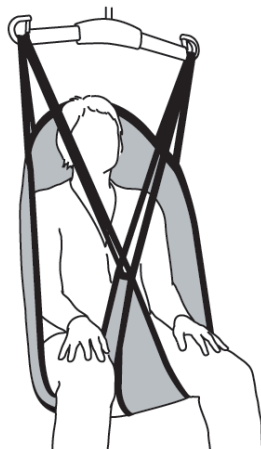ShellHBSling is easy for the caregiver to apply and it gives the patient a feeling of security. It features leg sections, which means that the patient does not sit on the sling. Instead, the sling provides support from the tailbone up over the head and under the legs. ShellHBSling has a concave shape and provides excellent support for the head. The sling enables an upright seated posture that can be regulated via the two sling loop alternatives.

SystemRoMedic includes a wide range of functional, comfortable, high-quality slings that can be adapted for different types of lifting and for patients with different needs. The choice of model and material depends on the patient's needs and the transfer situation.

Functional inspection

Visual and mechanical inspection

Check the condition and function of the sling regularly. Always inspect the product after laundering. Check to ensure that seams and material are free from damage. Check to ensure that fabric is not worn or faded. Apply load to the device and check to ensure that clasps, handles, etc. withstand heavy load. If there are signs of wear, the product must be discarded.

Always read the manual

Always read the manuals for all assistive devices used during a transfer.

Keep the manual where it is accessible to users of the product.

Do not leave the patient unattended during a transfer situation.

To prevent discomfort and the risk of the patient sliding out of the sling, trial fit the sling carefully. First, apply the sling's back loops to the slingbar, and then apply the leg support loops. Use a slingbar that has been tested.

The lowest allowable safe working load always determines the safe working load of the assembled system.


Always check the safe working loads for the lift and accessories before use.

Contact your dealer if you have any questions.

Application of the sling when the patient is sitting in a wheelchair/chair or in bed

1

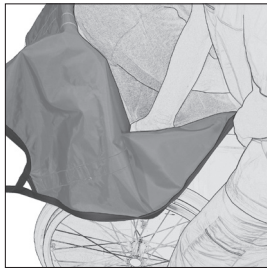2


1. Hold the sling in one hand and place your other hand in the pocket at the lower edge of the sling. This is easier if the patient leans forward or can be helped to lean forward slightly. Guide the sling down behind the patient's back, by sliding it along the back of the chair/wheelchair or along the bed, ideally, to the tailbone. The mid-section of the sling should follow the patient's spine.

2. Take care to guide the sling down to the patient's hips. Do this by pushing the sling into the space in the back corners of the wheelchair seat.

3. Then, work the leg sections out along the outside of the patient's thighs. Face the patient and pull the lower edge of the leg support, so that the sling slides into place around the hips. Light pressure against the patient's knees prevents him/her from sliding forward in the chair.


3

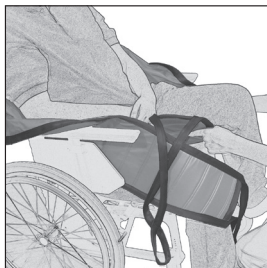4

4. After checking to ensure that the leg supports are of equal length on both sides, place them under the patient's thighs. Cross the leg supports and hook them to the slingbar. Start raising the lift, when the straps become taut, stop and check that everything is ok before proceeding with the lift.

Removal of the sling when the patient is sitting

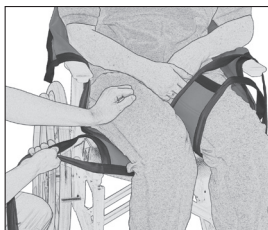1

2


1. Fold the leg support and allow it to slide under itself when you pull it out at the outside of the thigh.

2. Then, pull the sling up behind the patient's back.

Application of the sling when the patient is lying in bed, on a stretcher or on the floor


1

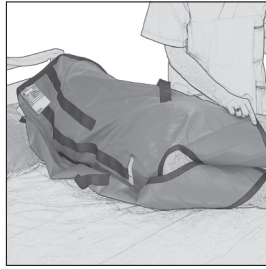2

1. Turn the patient towards you. Place the sling so that the product label is facing the underlying surface and the sling's mid-section follows the spine. The upper edge of the sling should be high enough so that the sling supports the head and the lower edge is level with the tailbone. Place the sling as far in under the patient as possible. Push the back-section loop under, so that you can reach it once you have turned the patient back onto the sling.

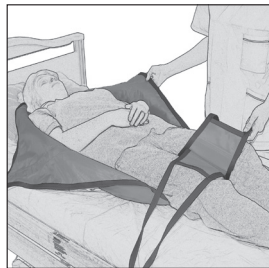3


4

2. To reduce the number of application steps, you can fold the back section up against the patient's back and place the leg supports between the patient's knees when he/she is turned towards you, and then pull them forward when the patient is lying on the sling.


3. Turn the patient back onto the sling and pull out the head support, back section and leg support.

4. Cross the leg supports and hook them to the slingbar. Start raising the lift, when the straps become taut, stop and check that everything is ok before proceeding with the lift.

Removal of the sling when the patient is recumbent

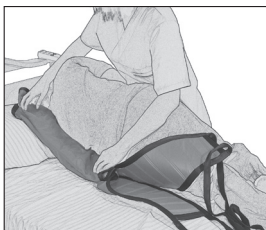1

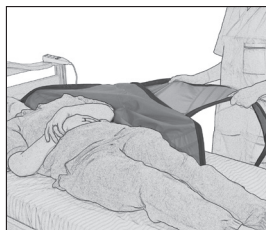
2

1. Lay the leg supports between the patient's knees.

2. Turn the patient towards you and fold the sling in under the patient so that the fabric slides towards you when you pull it out after you have turned back the patient.

Alternative application of the sling's leg supports

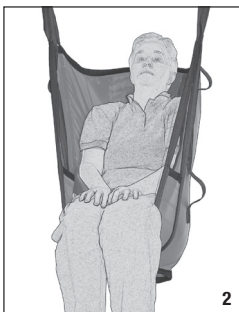1. In the standard application the leg supports are crossed before they are hooked to the slingbar.
 2. If it is desired to keep the legs of the patient together, for dignity or comfort the leg supports can also be overlapped under both legs.
 3. For more access to the perineal area of the patient. You can even choose not to cross the leg supports.
- Note, however, that this entails a risk that the patient might slide out of the sling.

Lifetime information & Periodic inspection

The product has an expected lifetime depending on how the product is used according to the manual instructions. The service life of the product varies depending on usage frequency, washing, drying, etc. Please be sure to read the product label and follow the care instructions for this product on www.directhealthcaregroup.com

The periodic inspection of the non-rigid body support unit should be performed at the timestamps stated by the manufacturer, but at least every 6 months to be sure. More frequently inspections may be required if a non-rigid body support unit is used or cleaned more frequently than advised by the manufacturer. The periodic inspection instructions for this product are found on www.directhealthcaregroup.com

Material:

Polyester: Durable material with low friction. Withstands high laundry temperatures.

Care of the product:

Read the product label.

Do not use rinsing agent.
To prolong product life,
avoid tumble-drying.

Size information


Slings can often be used for both 2 and 4-point suspension. It is essential to ensure that the width of the slingbar is appropriate for the sling size and the patient's status and capacity.

Use patient-specific slings to ensure safety and hygiene.

Size		XXS/Grey	XS/Orange	S/Red	M/Yellow	L/Green	XL/Purple	XXL/White
Rec. user weight	lbs	22-66	44-77	66-121	110-198	176-265	243-353	331-463
	kg	10-30	20-35	30-55	50-90	80-120	110-160	150-210
Sling overall height	inch	31,5	33,5	37	38	38,5	39,5	39,5
	cm	80	85	94	97	98	100	100
Sling seat measurement	inch	6,5	10	14	16	17,5	21,5	25,5
	cm	17	25	35	40	45	55	65

Contact your local distributor if you have any questions about the product and its use.

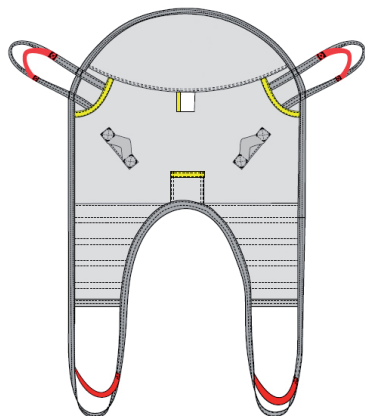See www.directhealthcaregroup.com for a complete list of distributors. Always make sure that you have the right version of the manual. The most recent editions of manuals are available for downloading from our website, www.directhealthcaregroup.com.

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Direct Healthcare Group Sverige AB