Ready for life

## Manual - Raiser bed and chair

Manual nr: 00761- EN, Nord Ver. 1 110330

Ready for life is a series of assistive devices that are designed to give users the best possible conditions for living an active life on a daily basis. The products have been specially developed to meet the conditions for homecare and institutional care settings. Ready for life offers the right quality, the right function and the right performance at the right price.

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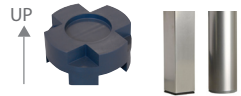
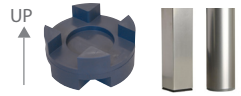
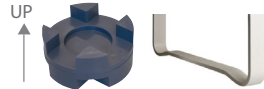
### English

By using 1-3 furniture raisers under each corner of the bed, the bed can be raised up to 150 mm depending on what side of the raiser is turned to face upwards.

### For beds with round or square legs

One end of the raiser has a raised cross with an 8 cm Ø round recess in the middle. This side is turned to face upwards for use with beds with round or square legs. The bed legs are placed in the recess in the middle, which raises the bed by 60 mm. By stacking several raisers on top of each other, the bed can be raised an additional 45 mm for each set of raisers.

### For beds with round or square legs, or beds with runners

The other end of the raiser has a recessed cross with an 8 cm Ø round recess in the middle. This end is turned to face upwards to enable raising of beds with round or square legs, or beds with runners. The bed legs are placed in the recess in the middle, which raises the bed by 35 mm. Runners are placed in the re-cessed cross, which has two different widths to enable optimal fit and stability. Beds with runners can be raised at 45 mm intervals. By stacking several raisers on top of each other, the bed can be raised an additional 45 mm for each set of raisers.

The furniture raisers can also be used to tilt the bed by raising the height at one end. This is done by installing raisers on one end of the bed only. Never use more than one set of raisers to tilt a bed.

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## English

Ready for life also features a double-sided furniture raiser that can be used to raise the height of beds and armchairs to a range of height positions. Which end of the raiser faces up or down, respectively, depends on whether the bed frame has legs or runners and the height to which the bed is to be raised.

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## Control of function

### Visual inspection

Test the function of the product regularly. Check to ensure that there is no material damage.

### Before use:

Check to ensure that the product is installed correctly.

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## Always read the manual

Always read the manual.  
Keep the manual where it is easily accessible for users of the product.

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### English

### Important information

• Raising furniture always affects the stability. An individual assessment regarding the center of mass and the furniture's support face as well as the movement patterns of the user must always be done.

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## Maintenance

Clean Raiser bed and chair with a mild, non-abrasive, non-corrosive cleaning agent. Do not use cleaning agents that contain phenol or chlorine, as this may damage the material.

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