
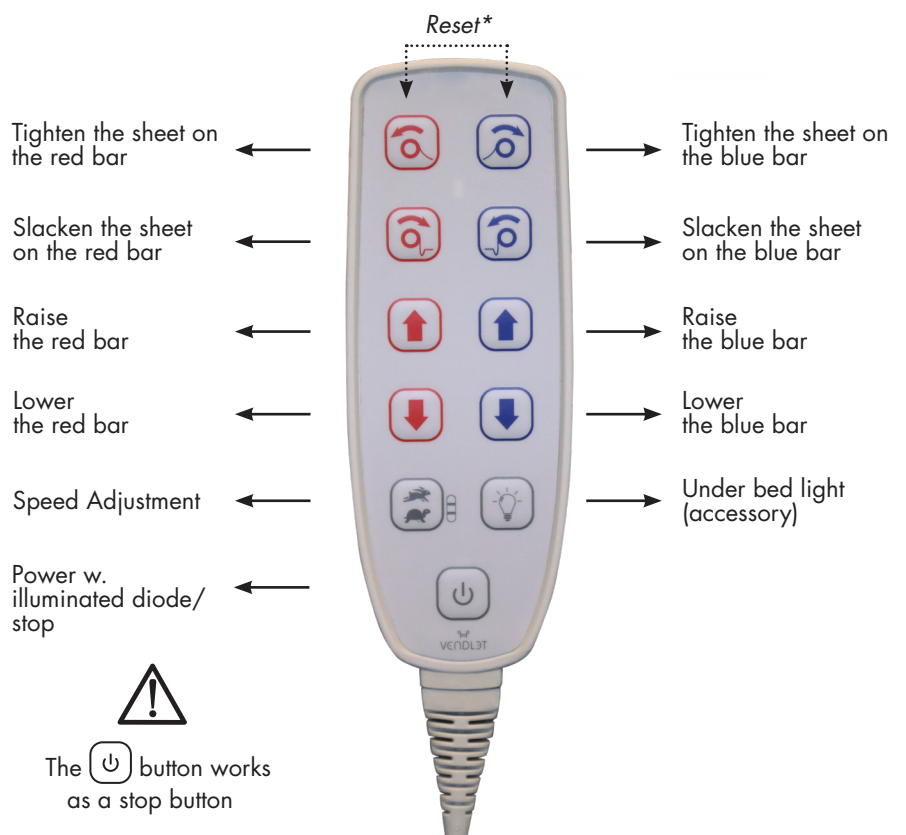


VENDLET V5S Speed Adjust VENDLET Bari

QUICK GUIDE

THE USE

1. Elevate the bed to the right working height.
2. Double-click the  button to activate the VENDLET system. The VENDLET system produces an acoustic signal and the green diode over the power symbol turns on. When the hand control has been inactive for 90 seconds, the system is off again.
3. The patient is repositioned by operating the hand control.



4. If it is estimated that the bedridden must use side rails, raise the bars. Slacken the sheet.
5. If it is estimated that the bedridden should not use side rails, lower the bars. Slacken the sheet.

Always, keep focus on the client when using VENDLET.



* Reset the system by simultaneously pressing the two top buttons on the hand control. You will hear 10 rhythmic beeps. Press and hold the buttons until the beeping stops.

**VENDLET must be used
and installed in accordance
with the manual**

Scan the QR code
to watch instruction
videos



VENDLET V5S Speed Adjust VENDLET Bari

QUICK GUIDE



ATTACHING THE SLIDE SHEET

1. Fold the slide sheet corner elastic straps around the mattress.
2. Check that the surface is smooth and even without creases and folds. Make sure that slippery side is facing upwards.
3. Fasten the velcro tapes around the mattress frame.
It is important that the tapes are fastened around the frame under the mattress and not around the bed frame as this will prevent the back, thigh, and foot parts of the bed from being raised.







ATTACHING THE TURNING SHEET

1. Place the turning sheet in the middle of the bed, velcro tape facing downward.
2. Press the turning sheet onto the velcro strip on one of the bars; begin at the head of the bed and work your way down to the foot of the bed.
3. The turning sheet must always be placed OVER the bars to prevent the fingers, hair ect. of the patient from getting caught under the sheet.



The pulling power will be reduced and the motors will switch off when the turning sheet is pulled in by the bar if the sheet is placed under the bars.

4. Press the  or  button on the hand control and allow approx. two-thirds of the sheet to be pulled in by the bar.
5. Press the turning sheet onto the velcro strip on the opposite bar.
6. Use this bar to pull in the turning sheet by pressing the  or  button until the sheet covers the mattress evenly with the red line at the centre of the mattress.



*VENDLET must be used
and installed in accordance
with the manual*

Scan the QR code
to watch instruction
videos

