

XXL X'PRESS BED



THIS PRODUCT MEETS THE REQUIREMENTS OF EUROPEAN DIRECTIVE 93/42/CEE IT HAS BEEN DESIGNED AND MADE IN ACCORDANCE WITH STANDARD CEI60601-2-52



The first extra-wide free-standing bed without a transport kit.



• Available in **120 and 140 cm**

• **Quick assembly** in minutes, without tools




• **Extremely manoeuvrable:** the bed can be moved around using its 100 mm diameter wheels, which optimise its ability to clear obstacles.




Ideal for at home
The XXL X'PRESS bed can be put upright without tools



The XXL X'PRESS bed is available in 2 versions: 2 and 3 electric functions with folding.

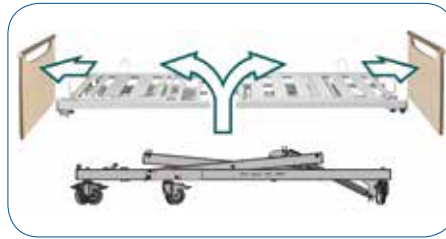

Remove the panelling, the backrest and then the leg-rest.

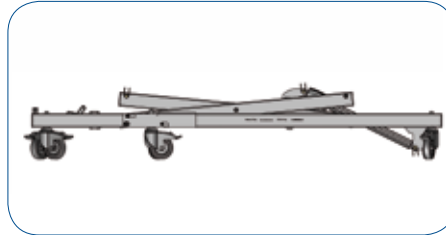


Slightly raise the crosspiece to raise the base.

THE XXL X'PRESS BED: CONFIGURABLE AND ADAPTABLE (available options)



Electric backrest (original).

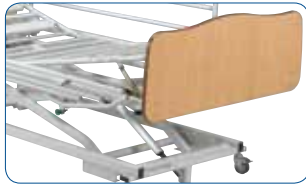

Manual leg-rest using racks.



Wheels Ø 100 mm.

Independent breaking for each bed wheel.



Large welded metal slats. These are designed for easy maintenance.



Base extension of 20 cm

TECHNICAL CHARACTERISTICS

PRODUCT NAME	SIZES	RB	LEG-REST	STANDARD USE LOADING	WHEELS	STANDARDS	WATERPROOF
XXL X'PRESS 2 electric functions	120 x 200 cm	70°	12° racks	315 kg 30 kg mattress 270 kg patient 15 kg accessories	Ø 100 mm with breaks without wire guard	CEI 60601-2-52	LINAK IP66 remote, box, motor V0
	140 x 200 cm						
XXL X'PRESS 3 electric functions	120 x 200 cm		Knee bend of 26° (thighs and 16° foot adjustment using racks)				
	140 x 200 cm						
Accessories	Epoxy offset trapeze bar (see accessory sheet)						
	VISION epoxy folding rails (see accessory sheet)						
	Abélia head/baseboards, 120 or 140 cm (see panel and colour chart sheet)						

- 2 trapeze supports (at the head of the bed).
- 2 accessory supports (on the long sections).
- 2 drip stand supports (at the head of the bed).

5 year* mechanical, electrical and remote guarantee. 1 year* panel and accessory guarantee. *Under standard conditions of use and excluding wear and tear - see instructions for use.

MÉDICATLANTIC Groupe Winncare - 85670 Saint-Paul-Mont-Penit - France
Tel : +33 (0)2 51 98 55 64 - Fax: +33 (0)2 51 98 59 07 - E-mail : info@medicatlantic.com

www.winncare.fr

With a focus on the continuous improvement of its products, MEDICATLANTIC reserves the right to modify technical characteristics without prior notice - non-contractual photos.

FP000010



WINN CARE
GROUP

by

MÉDICATLANTIC