

X'PRESS BED

2nd generation

THIS PRODUCT SATISFIES THE REQUIREMENTS OF EUROPEAN DIRECTIVE 93/42/CEE. IT IS DESIGNED AND MANUFACTURED IN COMPLIANCE WITH THE CEI 60601-2-52 STANDARDS IN FORCE.

• **Ergonomics:** The mechanical ergonomics of the bed enable the installer to adopt the positions recommended for preventing MSD (musculoskeletal disorders) during the delivery, assembly and disassembly of the bed.

Its small overall dimensions and storage on its own wheels facilitate transfers in all situations (lifts, small vehicles, etc).

This self-supporting bed does not require a trolley.

• **Simplicity:** The assembly and disassembly operations are very implicit.

No tools are required.

Blue coloured parts are associated with functions. Their role is to guide bed assembly operations.

• **Safety:** The new "Easy Move" panel mounting system prevents the loss of rondos and makes the panel mountings on the bed more secure.

Cable routing is designed to reduce tension and the risk of damage.



Self-supporting 2nd generation X'PRESS BED



Conversion of self-supporting trolley to become the chassis




Locking of the 1/2 base with the legs horizontal.



Fixing together of the two 1/2 bases by means of through-tube rondos

EASY MOVE

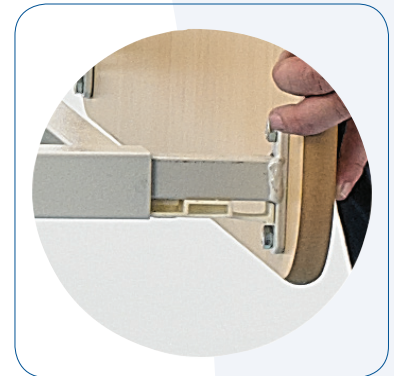

Panel mounting by the EASY MOVE system: simple, fast and secure.

TECHNICAL CHARACTERISTICS

DESCRIPTION	DIMENSIONS	VARIABLE HEIGHT	BACK SUPPORT	LEG SUPPORT	RATED WORKING LOAD	WHEELS	STANDARDS	SEALING
Two-function 2 nd generation X'PRESS	Assembled bed: 90 x 200 cm, dividable base. Stored bed: 50 cm wide, 90 cm long and 160 cm high.	34 to 84 cm	0 to 75°	12° ratchet	170kg Mattress 20 kg Patient 135 kg Accessories 15 kg	100 mm dia. with brakes and without wire protector	CEI 60601-2-52	Motor and box, LINAK IP66 Lockable remote control, LINAK IP54
Three-function 2 nd generation X'PRESS without knee bending				13.5° electric				
Three-function 2 nd generation X'PRESS with knee bending				Knee bending of 26° and 16° of foot lifting by ratchet.				

- 2 accessory brackets (on the long sides)
- 2 intravenous post support brackets (at the bed head)
- 2 bed gallows brackets (at the bed head)

Warranty:
Mechanical structure, jacks and remote control - 5 years*
Panels and accessories - 1 year*
* In normal usage conditions and excluding wearing parts – see instructions in user's manual

MÉDICATLANTIC Groupe Winncare - 85670 Saint-Paul-Mont-Penit - France
Phone: +33 (0)2 51 98 55 64 - Fax: +33 (0)2 51 98 59 07 - E-mail: info@medicatlantinc.fr

www.winncare.fr



Within their policy of constant product improvement
MÉDICATLANTIC reserve the right to make modifications to technical characteristics without prior notice – Non-contractual photographs.

