

# Raizer® II

Lifting Chair

80 lifts  
per charge

Raizer II is a simple battery-operated mobile lifting chair that helps a fallen person up to an almost standing position in just 30 seconds. Raizer can easily be assembled and operated by only one person and does not require any physical effort besides a supportive hand. The design of the Raizer makes it easy to assemble and operate, and ensures easy cleaning.

## Key features

4 IDENTICAL LEGS and 2 IDENTICAL BACKRESTS ensure correct assembly every time


All SURFACE MATERIALS have been specified to be CLEANED EASILY and DISINFECT

ADDITIONAL SEAT BELTS can be fitted for increased stability of the fallen person if this is required

INTUITIVE SOUND and LED INDICATORS for correct assembly

DISPOSABLE HYGIENE COVERS are also available (optional accessory)

The REMOTE CONTROL is stored on the seat unit

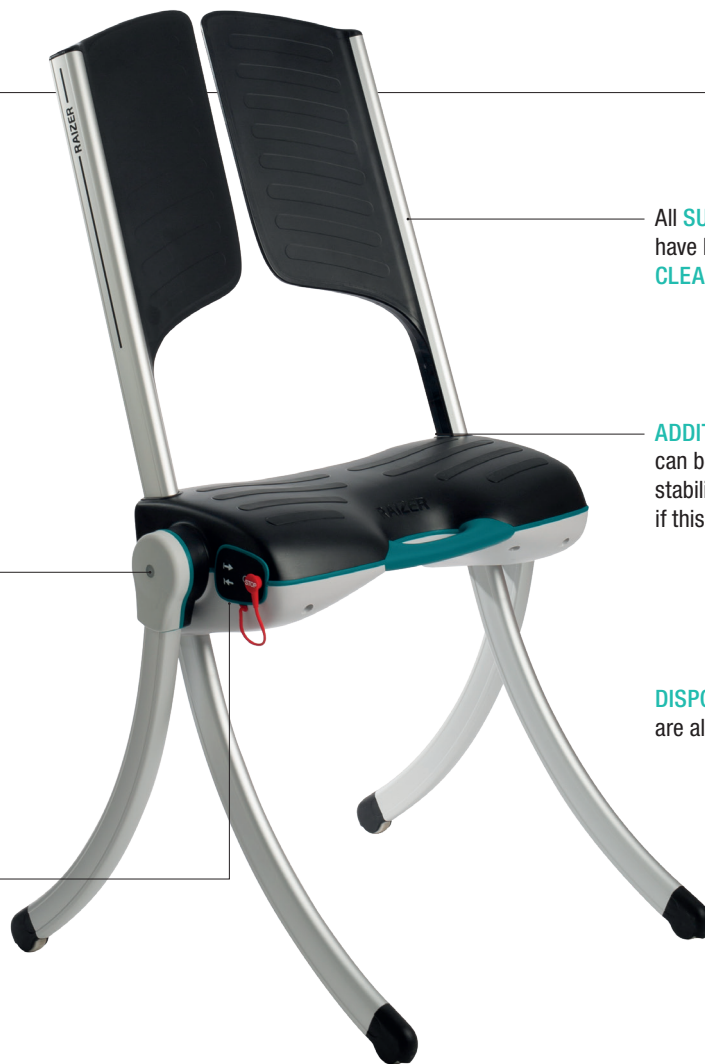


### Convenient and ergonomic to transport

- Effective transport
- All parts are supplied in two practical bags
- Legs, backrest and other parts are positioned intuitively in an easy-to-carry bag as their position in the bag mirrors the sequence of assembly
- The set unit is carried separately in a protective cover with a soft handle, ergonomically designed to make sure that it is carried close to the body
- The two bags can also be carried on a trolley (option accessory) with enough space to place the bags side by side

**DHG**  
Moving Health Forward

# Raizer® II

## Lifting Chair

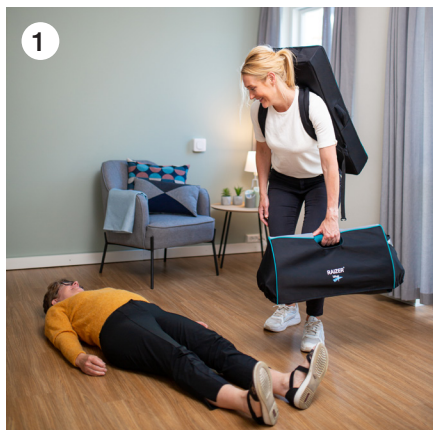

### Technical specification

Model	Raizer II
Lifting capacity	Max. 150kg
Lifting time	20-30 seconds
Weight	Seat: 8.5kg Legs and back support: 4.5kg Total: 13kg
Number of lifts when fully charged	40 at max. load. Approx. 80 with an average load
Charging time	Up to 6 hours
Charging of empty battery	10-15 min = 1 lift
Chargers	USB charger and 15W charger for standard power socket