

SystemRoMedic™

The essentials sling selection

A durable, functional range offering the user
a quality fit and comfortable experience



DHG
Moving Health Forward

quality functionality durability

SystemRoMedic™ Slings are characterised by a quality fit and functionality, as well as practical, but not unnecessary features.

All details are placed where they provide the most benefit, without being in the way for either the user or the caregiver.


- Durable fabrics with maximum service life:
- Strong seams, reinforced in strategic places
- Strong edging for better wear
- All belts and fasteners are strong, safe and easy to adjust
- Fabrics that are easy to care for and withstand laundering at high temperatures
- Disposable material for short-term, patient specific use

NEW PROMO PRICE Contact your local DHG representative for a quote today







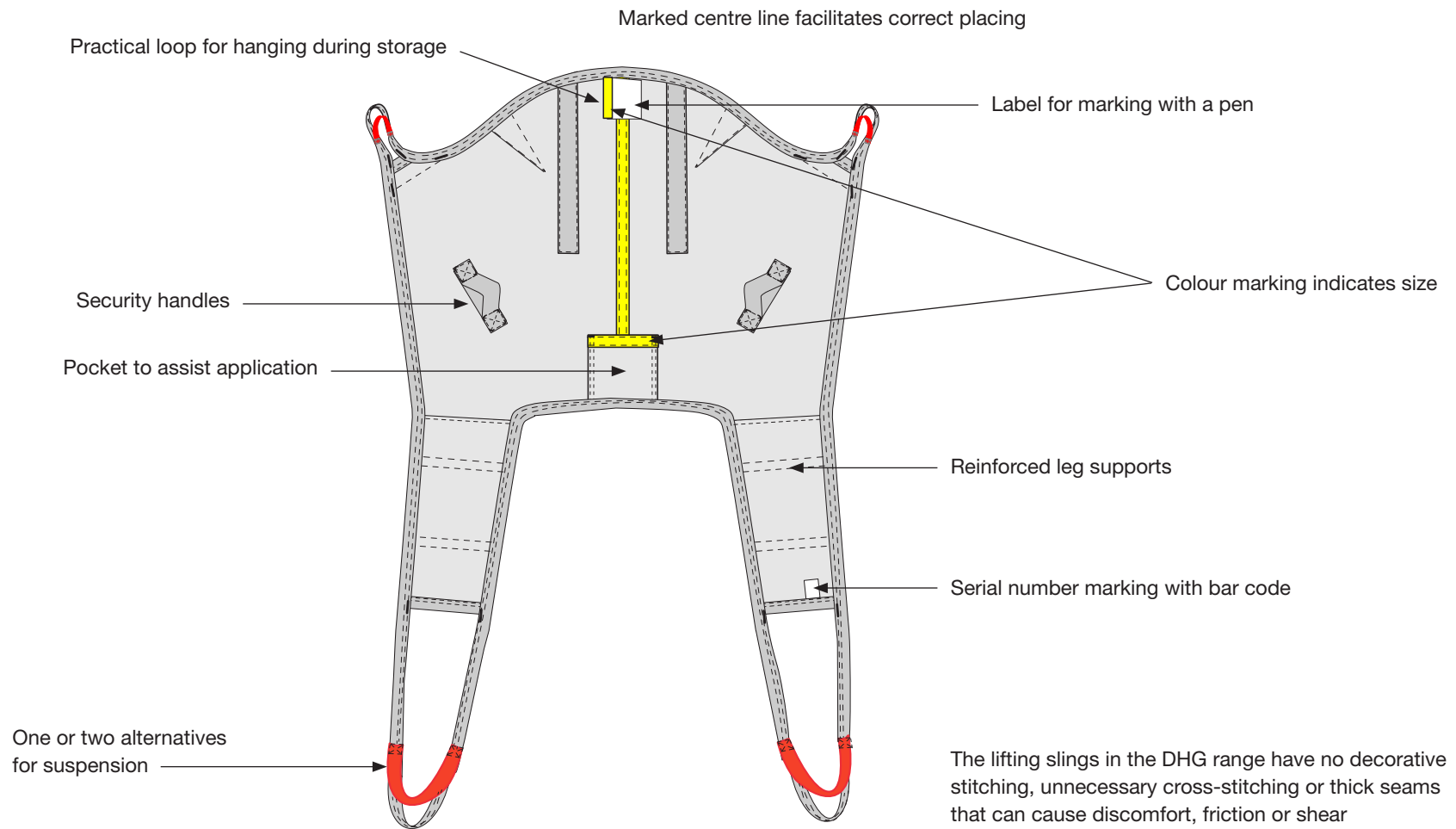








Simplicity is comfort




Different models for different situations

DHG offers solutions for all different transfer situations and for all users. Each of our lifting sling models is designed and produced for optimal functionality, both for the user and the caregiver, for one or a few specific situations. Some models and materials work best for seated transfers, for example, between bed and chair, while others are adapted for situations such as lifting from the recumbent position, for hygiene, bathing or showering. The table to the right, shows which models are most suitable in different situations.

		Situations										
												
		Supine transfer	Turning in bed	Higher up in bed	Up from the floor	Support for gait training	Standing training	From bed to wheelchair	From bed to commode	From wheelchair to toilet	Rising – sitting down	Positioning in chair
Seated lifting, divided leg supports	Lifting Slings											
	HighBackSling				■			■				■
	ShellHBSling				■			■				■
	FlexibleSling				■			■				■
	HygieneHBSling								■	■		■
	AmpHBSling				■			■				■
	BasicSling							■				■
	AmpSling							■				■
	LowBackSling							■				■
	HygieneSling								■	■		
Seated lifting, undivided leg supports	HygieneLBSling								■	■		■
	FlexibleSling, undivided							■				■
	ClassicHBSling							■				■
	ClassicSling							■				■
	ClassicSling, commode							■	■	■		■
	FlexibleSling, hygiene							■	■	■		■
	BariSling							■				
Horizontal lifting/Turning	StretcherSling	■		■	■							
	SafeHandlingSheet	■	■	■	■							
Rising/ Standing/ Walking	ThoraxSling						■	■	■	■	■	■
	StandingVest						■				■	
	WalkingVest					■	■				■	
	BariVest					■	■			■		

Size is everything

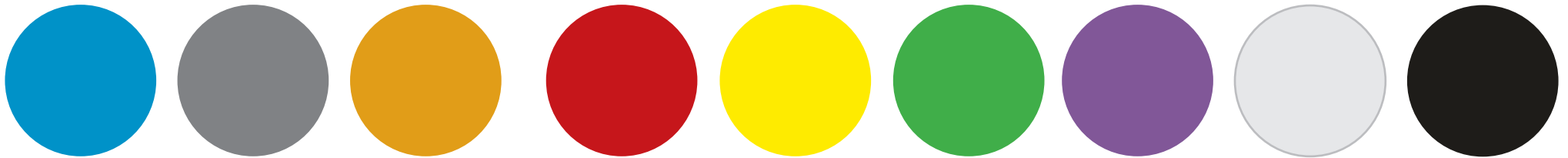

DHG's lifting slings are available in sizes ranging from XXS to XXL. The smallest sizes fit children and small adults, and the largest sizes fit large and, often, overweight users. The three sizes in the middle of the scale, S, M and L, fit most normal-sized adults. We also have special lifting sling models in XXXL and XXXXL for extremely large and heavy users.

Trial fitting is a must

The size guide measurements and recommendations are only intended as an aid to selecting lifting slings.

For a transfer with a stationary or mobile lift to be as safe and comfortable as possible, it is very important that the lifting sling is the right size and that the model is appropriate for the particular user and lifting situation.


To determine which model, size and material are most suitable for the individual user, individual trial fitting is always necessary.

The lifting sling is too large:

- if the user's seat slides out too much in the lower opening
- if the leg supports slide forward towards the hollows of the user's knees

The lifting sling is too small:

- if the user's face is too close to the sling bar
- if the leg supports chafe at the groin
- if the lifting sling looks or feels too tight

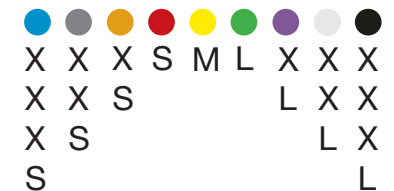

Well-established colour coding system

We use an international colour coding system, with different colours for different sizes, to facilitate choice of size.

The colour marking is affixed to the lifting sling in the form of a coloured ribbon that also functions as a loop that can be used to hang up the sling. This allows you to find the right size at a glance.

Durable materials



Durable, easy-to-care-for polyester

DHG's lifting slings in polyester are comfortable, durable and easy to maintain. Polyester is a strong and durable material. It is also very pliable and has low friction, which makes the lifting sling easy to apply and remove, as well as soft and comfortable for the user.


The lifting slings have a repellent surface, but when they need to be laundered, they withstand high temperatures, which contributes to good hygiene.



Flexible, airy polyester net

Polyester net makes the lifting slings flexible, pliable and airy, which gives pain-sensitive users and users with special needs better comfort without jeopardizing safety or practical features.

DHG's polyester net is soft and pliable, allowing both air and moisture to pass through the mesh. It is a tough, easy-to-maintain material that withstands high laundry temperatures and can be left under the user between transfers.

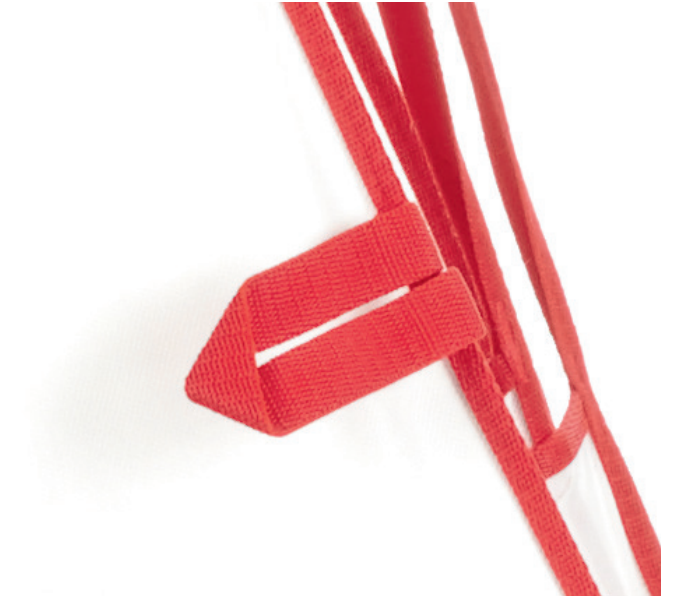

Disposable material for optimal hygiene

DHG also has lifting slings for patient specific use. These are made of a sturdy, durable non-woven material. Instead of laundering, these lifting slings are discarded after use, when they become soiled or when the user no longer needs them. This prevents the spread of multi-resistant and other types of infectious bacteria. Disposable material is a very practical solution which also enables optimal hygiene and safety.



BasicSling



Standard model for users with good head stability

BasicSling is designed for use in most common lifting situations and when the patient has poor torso stability but good head stability.

- Divided leg supports; easy for the caregiver to apply and remove, the user does not sit on the lifting sling
- Shoulder-high back; the user is supported from the tailbone up to the shoulders and under the legs
- Upright seated posture
- Standard model, suits most users with good head stability
- BasicSling in polyester net is available with or without reinforced leg supports
- Pocket at lower edge for easier application
- SWL 300kg/660lbs

HighBackSling



Standard model for users who need a lot of support

HighBackSling is designed for use in most common lifting situations and when the user needs extra head support. HighBackSling can also be used for lifting from a recumbent position on the floor.

- Divided leg supports; easy for the caregiver to apply and remove, the user does not sit on the lifting sling
- High back; the user is supported from the tailbone up over the head and under the legs
- Slightly reclined seated posture
- Standard model, suits most users
- Removable plastic stays give the head support extra stability
- HighBackSling in polyester net is available with or without reinforced leg supports
- Pocket at lower edge for easier application
- SWL 300kg/660lbs (disposable material: SWL 200kg/440lbs)

FlexibleSling



Pliable support for users with differing postural needs

FlexibleSling is a sling that meets the need for a lot of support and, in a flexible way, follows the contours of the body and gives the user a sense of security. FlexibleSling is designed for use in most common lifting situations, especially for users with postural needs and users who need extra head support. FlexibleSling can also be used for lifting from a recumbent position on the floor.

FlexibleSling is also available with undivided leg supports and with undivided leg supports and commode opening.

- Divided leg supports; easy for the caregiver to apply and remove, the user does not sit on the lifting sling
- High back; the user is supported from the tailbone up over the head and under the legs
- Slightly reclined seated posture
- Made of airy polyester net and can therefore be left in a wheelchair, moulded seat or similar after completed transfer
- Follows the shape of the body
- Particularly suitable for users in moulded seats
- FlexibleSling with divided leg supports is available with or without reinforced leg supports
- SWL 300kg/660lbs

HygieneSling










Three versions for more active users

HygieneSling is a lifting sling that is easy to apply and makes it easier for the caregiver to dress/undress the user in connection with toileting. HygieneSling, which is designed for users with adequate muscular tonus, is available in three versions; with or without a belt, or with a belt and extra back support.

- Divided leg supports; easy for the caregiver to apply and re-move, the user does not sit on the lifting sling
- Low back; the user is supported behind the back, up to the armpits and under the legs
- Slightly reclined seated posture
- The user must have adequate muscular tonus
- HygieneSling is available in three versions; with or without a belt, and with a belt and extra back support
- SWL 300kg/660lbs

Size guide for lifting slings



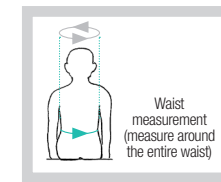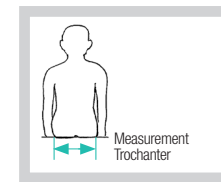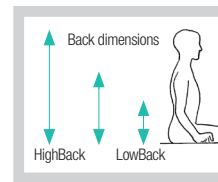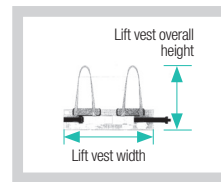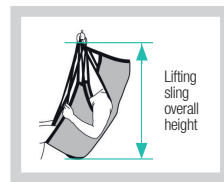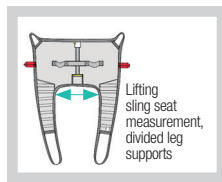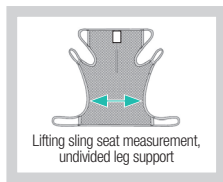



- To determine the most appropriate model, size and material for the individual user, trial fitting of lifting slings is always essential.
- The size guide's measurements and recommendations are only intended as an aid to selecting lifting slings.
- The lifting sling's seat width gives an indication of how wide the user's hips (trochanter width) can be. The user's trochanter width must be within +/- 5cm/2" of the lifting sling's seat width.
- The lifting sling is too large if the user's seat slides out too far at the lower opening, or if the leg supports slide forward towards the hollows of the user's knees.
- The lifting sling is too small if the user's face is too close to the slingbar, if the leg supports chafe at the groin, or if the sling appears to be too tight or feels too tight.



Divided leg supports/High back

	Art. No.	Product name	Material/Description	SWL	Size	Rec user weight	Sling seat measurement	Sling overall height
	4080002	HighBackSling	polyester	300kg/660lbs	XXS	10–25kg/22–55lbs	20cm/8"	77cm/30"
	4080003	HighBackSling	polyester	300kg/660lbs	XS	20–35kg/44–77lbs	25cm/10"	85cm/33"
	4080004	HighBackSling	polyester	300kg/660lbs	S	30–55kg/66–122lbs	35cm/14"	91cm/36"
	4080006	HighBackSling	polyester	300kg/660lbs	M	50–90kg/110–198lbs	45cm/18"	94cm/37"
	4080007	HighBackSling	polyester	300kg/660lbs	L	80–110kg/176–242lbs	56cm/22"	98cm/39"
	4080008	HighBackSling	polyester	300kg/660lbs	XL	110–160kg/242–353lbs	62cm/24"	98cm/39"
	4080009	HighBackSling	polyester	300kg/660lbs	XXL	130–200kg/287–441lbs	70cm/28"	102cm/40"
	4082002	HighBackSling	polyester net	300kg/660lbs	XXS	10–25kg/22–55lbs	20cm/8"	77cm/30"
	4082003	HighBackSling	polyester net	300kg/660lbs	XS	20–35kg/44–77lbs	25cm/10"	85cm/33"
4082004	HighBackSling	polyester net	300kg/660lbs	S	30–55kg/66–122lbs	35cm/14"	91cm/36"	
4082006	HighBackSling	polyester net	300kg/660lbs	M	50–90kg/110–198lbs	45cm/18"	94cm/37"	
4082007	HighBackSling	polyester net	300kg/660lbs	L	80–110kg/176–242lbs	56cm/22"	98cm/39"	
4082008	HighBackSling	polyester net	300kg/660lbs	XL	110–160kg/242–353lbs	62cm/24"	98cm/39"	
4082009	HighBackSling	polyester net	300kg/660lbs	XXL	130–200kg/287–441lbs	70cm/28"	102cm/40"	



Divided leg supports/High back (continued from previous page)

Art. No.	Product name	Material/Description	SWL	Size	Rec user weight	Sling seat measurement	Sling overall height	
	40822002	HighBackSling	polyester net, <i>not reinforced</i>	300kg/660lbs	XXS	10–25kg/22–55lbs	20cm/8"	77cm/30"
	40822003	HighBackSling	polyester net, <i>not reinforced</i>	300kg/660lbs	XS	20–35kg/44–77lbs	25cm/10"	85cm/33"
	40822004	HighBackSling	polyester net, <i>not reinforced</i>	300kg/660lbs	S	30–55kg/66–122lbs	35cm/14"	91cm/36"
	40822006	HighBackSling	polyester net, <i>not reinforced</i>	300kg/660lbs	M	50–90kg/110–198lbs	45cm/18"	94cm/37"
	40822007	HighBackSling	polyester net, <i>not reinforced</i>	300kg/660lbs	L	80–110kg/176–243lbs	56cm/22"	98cm/39"
	40822008	HighBackSling	polyester net, <i>not reinforced</i>	300kg/660lbs	XL	110–160kg/242–353lbs	62cm/24"	98cm/39"
	40822009	HighBackSling	polyester net, <i>not reinforced</i>	300kg/660lbs	XXL	130–200kg/287–441lbs	70cm/28"	102cm/40"
	40890004	HighBackSling	disposable, <i>non-woven</i>	200kg/440lbs	S	30–55kg/66–122lbs	35cm/14"	91cm/36"
	40890006	HighBackSling	disposable, <i>non-woven</i>	200kg/440lbs	M	50–90kg/110–198lbs	45cm/18"	94cm/37"
	40890007	HighBackSling	disposable, <i>non-woven</i>	200kg/440lbs	L	80–110kg/176–242lbs	56cm/22"	98cm/39"
	40890008	HighBackSling	disposable, <i>non-woven</i>	200kg/440lbs	XL	110–160kg/242–353lbs	62cm/24"	98cm/39"
	40890009	HighBackSling	disposable, <i>non-woven</i>	300kg/660lbs	XXL	150–250kg/330–551lbs	70cm/28"	102cm/40"
	40700002	ShellHBSling	polyester	300kg/660lbs	XXS	10–30kg/22–66lbs	17cm/7"	80cm/31"
	40700003	ShellHBSling	polyester	300kg/660lbs	XS	20–35kg/44–77lbs	25cm/10"	85cm/33"
	40700004	ShellHBSling	polyester	300kg/660lbs	S	30–55kg/66–122lbs	35cm/14"	94cm/37"
	40700006	ShellHBSling	polyester	300kg/660lbs	M	50–90kg/110–198lbs	40cm/16"	97cm/38"
	40700007	ShellHBSling	polyester	300kg/660lbs	L	80–120kg/176–265lbs	45cm/18"	98cm/39"
	40700008	ShellHBSling	polyester	300kg/660lbs	XL	110–160kg/242–353lbs	55cm/22"	100cm/39"
	40700009	ShellHBSling	polyester	300kg/660lbs	XXL	150–210kg/331–463lbs	65cm/26"	100cm/39"
	42520002	FlexibleSling	polyester net	300kg/660lbs	XXS	10–30kg/22–66lbs	20cm/8"	75cm/30"
	42520003	FlexibleSling	polyester net	300kg/660lbs	XS	20–35kg/44–77lbs	25cm/10"	81cm/32"
	42520004	FlexibleSling	polyester net	300kg/660lbs	S	25–55kg/55–121lbs	35cm/14"	90cm/35"
	42520006	FlexibleSling	polyester net	300kg/660lbs	M	50–75kg/110–165lbs	45cm/18"	102cm/40"
	42520007	FlexibleSling	polyester net	300kg/660lbs	L	70–110kg/154–242lbs	55cm/22"	102cm/40"
	42520008	FlexibleSling	polyester net	300kg/660lbs	XL	100–130kg/220–287lbs	65cm/26"	104cm/41"
	42522002	FlexibleSling	polyester net, <i>not reinforced</i>	300kg/660lbs	XXS	10–30kg/22–66lbs	20cm/8"	75cm/30"
	42522003	FlexibleSling	polyester net, <i>not reinforced</i>	300kg/660lbs	XS	20–35kg/44–77lbs	25cm/10"	81cm/32"
	42522004	FlexibleSling	polyester net, <i>not reinforced</i>	300kg/660lbs	S	25–55kg/55–121lbs	35cm/14"	90cm/35"
	42522006	FlexibleSling	polyester net, <i>not reinforced</i>	300kg/660lbs	M	50–75kg/110–165lbs	45cm/18"	102cm/40"
	42522007	FlexibleSling	polyester net, <i>not reinforced</i>	300kg/660lbs	L	70–110kg/154–242lbs	55cm/22"	102cm/40"
	42522008	FlexibleSling	polyester net, <i>not reinforced</i>	300kg/660lbs	XL	100–130kg/220–287lbs	65cm/26"	104cm/41"
	42800004	AmpHBSling	polyester	300kg/660lbs	S	35–60kg/77–132lbs	37cm/15"	90cm/35"
	42800006	AmpHBSling	polyester	300kg/660lbs	M	50–80kg/110–176lbs	45cm/18"	97cm/38"
	42800007	AmpHBSling	polyester	300kg/660lbs	L	70–110kg/154–242lbs	50cm/20"	100cm/39"
	42800008	AmpHBSling	polyester	300kg/660lbs	XL	100–130kg/220–286lbs	60cm/24"	100cm/39"






Divided leg supports/Shoulder-high back

	Art. No.	Product name	Material/Description	SWL	Size	Rec user weight	Sling seat measurement	Sling overall height
	40500002	BasicSling	polyester	300kg/660lbs	XXS	10–25kg/22–55lbs	21cm/8"	63cm/25"
	40500003	BasicSling	polyester	300kg/660lbs	XS	20–35kg/44–77lbs	27cm/11"	65cm/26"
	40500004	BasicSling	polyester	300kg/660lbs	S	35–55kg/77–121lbs	37cm/15"	77cm/30"
	40500006	BasicSling	polyester	300kg/660lbs	M	45–80kg/99–176lbs	50cm/20"	89cm/35"
	40500007	BasicSling	polyester	300kg/660lbs	L	75–110kg/165–242lbs	60cm/24"	95cm/37"
	40500008	BasicSling	polyester	300kg/660lbs	XL	100–150kg/220–330lbs	65cm/26"	98cm/39"
	40500009	BasicSling	polyester	300kg/660lbs	XXL	140–200kg/309–441lbs	75cm/30"	98cm/39"
	40520002	BasicSling	polyester net	300kg/660lbs	XXS	10–25kg/22–55lbs	21cm/8"	63cm/25"
	40520003	BasicSling	polyester net	300kg/660lbs	XS	20–35kg/44–77lbs	27cm/11"	65cm/26"
	40520004	BasicSling	polyester net	300kg/660lbs	S	35–55kg/77–121lbs	37cm/15"	77cm/30"
	40520006	BasicSling	polyester net	300kg/660lbs	M	45–80kg/99–176lbs	50cm/20"	89cm/35"
	40520007	BasicSling	polyester net	300kg/660lbs	L	75–110kg/165–242lbs	60cm/24"	95cm/37"
	40520008	BasicSling	polyester net	300kg/660lbs	XL	100–150kg/220–330lbs	65cm/26"	98cm/39"
	40520009	BasicSling	polyester net	300kg/660lbs	XXL	140–200kg/309–441lbs	75cm/30"	98cm/39"
	40522002	BasicSling	polyester net, <i>not reinforced</i>	300kg/660lbs	XXS	10–25kg/22–55lbs	21cm/8"	63cm/25"
	40522003	BasicSling	polyester net, <i>not reinforced</i>	300kg/660lbs	XS	20–35kg/44–77lbs	27cm/11"	65cm/26"
	40522004	BasicSling	polyester net, <i>not reinforced</i>	300kg/660lbs	S	35–55kg/77–121lbs	37cm/15"	77cm/30"
	40522006	BasicSling	polyester net, <i>not reinforced</i>	300kg/660lbs	M	45–80kg/99–176lbs	50cm/20"	89cm/35"
	42700003	AmpSling	polyester	300kg/660lbs	XS	20–40kg/44–88lbs	27cm/11"	85cm/33"
	42700004	AmpSling	polyester	300kg/660lbs	S	35–60kg/77–133lbs	37cm/15"	90cm/35"
	42700006	AmpSling	polyester	300kg/660lbs	M	50–80kg/110–176lbs	45cm/18"	97cm/38"
	42700007	AmpSling	polyester	300kg/660lbs	L	70–110kg/154–242lbs	50cm/20"	100cm/39"
	42700008	AmpSling	polyester	300kg/660lbs	XL	100–130kg/220–287lbs	60cm/24"	100cm/39"
	42700009	AmpSling	polyester	300kg/660lbs	XXL	120–180kg/265–397lbs	70cm/28"	100cm/39"
	42720003	AmpSling	polyester net	300kg/660lbs	XS	20–40kg/44–88lbs	27cm/11"	85cm/33"
	42720004	AmpSling	polyester net	300kg/660lbs	S	35–60kg/77–133lbs	37cm/15"	90cm/35"
	42720006	AmpSling	polyester net	300kg/660lbs	M	50–80kg/110–176lbs	45cm/18"	97cm/38"
	42720007	AmpSling	polyester net	300kg/660lbs	L	70–110kg/154–242lbs	50cm/20"	100cm/39"
	42720008	AmpSling	polyester net	300kg/660lbs	XL	100–130kg/220–287lbs	60cm/24"	100cm/39"
	42720009	AmpSling	polyester net	300kg/660lbs	XXL	120–180kg/265–397lbs	70cm/28"	100cm/39"

Divided leg supports/Low back

	Art. No.	Product name	Material	SWL	Size	Rec user weight	Sling seat measurement	Sling overall height	
	4020002	LowBackSling	polyester	300kg/660lbs	XXS	10–25kg/22–55lbs	21cm/8"	60 cm/24"	
	4020003	LowBackSling	polyester	300kg/660lbs	XS	20–35kg/44–77lbs	27cm/11"	60 cm/24"	
	4020004	LowBackSling	polyester	300kg/660lbs	S	30–50kg/66–110lbs	37cm/15"	77 cm/30"	
	4020006	LowBackSling	polyester	300kg/660lbs	M	45–80kg/99–176lbs	47cm/18"	85 cm/33"	
	4020007	LowBackSling	polyester	300kg/660lbs	L	75–110kg/165–242lbs	57cm/22"	90 cm/35"	
	4020008	LowBackSling	polyester	300kg/660lbs	XL	100–140kg/220–309lbs	67cm/26"	92 cm/36"	
	4020009	LowBackSling	polyester	300kg/660lbs	XXL	130–190kg/287–419lbs	77cm/30"	92 cm/36"	
		4022002	LowBackSling	polyester net	300kg/660lbs	XXS	10–25kg/22–55lbs	21cm/8"	60 cm/24"
		4022003	LowBackSling	polyester net	300kg/660lbs	XS	20–35kg/44–77lbs	27cm/11"	60 cm/24"
4022004		LowBackSling	polyester net	300kg/660lbs	S	30–50kg/66–110lbs	37cm/15"	77 cm/30"	
4022006		LowBackSling	polyester net	300kg/660lbs	M	45–80kg/99–176lbs	47cm/18"	85 cm/33"	
4022007		LowBackSling	polyester net	300kg/660lbs	L	75–110kg/165–242lbs	57cm/22"	90 cm/35"	
4022008		LowBackSling	polyester net	300kg/660lbs	XL	100–140kg/220–309lbs	67cm/26"	92 cm/36"	
4022009		LowBackSling	polyester net	300kg/660lbs	XXL	130–190kg/287–419lbs	77cm/30"	92 cm/36"	

Undivided leg supports/High and shoulder-high back

	Art. No.	Product name	Material	SWL	Size	Rec user weight	Sling seat measurement	Sling overall height
	42521003	FlexibleSling, undivided	polyester net	300kg/660lbs	XS	20–40kg/44–88lbs	75cm/30"	80cm/31"
	42521004	FlexibleSling, undivided	polyester net	300kg/660lbs	S	30–70kg/66–154lbs	85cm/33"	85cm/33"
	42521006	FlexibleSling, undivided	polyester net	300kg/660lbs	M	60–90kg/132–198lbs	94cm/37"	93cm/37"
	42521007	FlexibleSling, undivided	polyester net	300kg/660lbs	L	80–130kg/176–287lbs	105cm/41"	96cm/38"
	42521008	FlexibleSling, undivided	polyester net	300kg/660lbs	XL	110–180kg/242–397lbs	115cm/45"	98cm/39"
		41522003	ClassicHBSling	polyester net	300kg/660lbs	XS	20–40kg/44–88lbs	30cm/12"
41522004		ClassicHBSling	polyester net	300kg/660lbs	S	30–60kg/66–132lbs	40cm/16"	98cm/39"
41522006		ClassicHBSling	polyester net	300kg/660lbs	M	55–85kg/121–187lbs	48cm/19"	106cm/42"
41522007		ClassicHBSling	polyester net	300kg/660lbs	L	80–120kg/176–265lbs	65cm/26"	110cm/43"
41522008		ClassicHBSling	polyester net	300kg/660lbs	XL	110–140kg/242–309lbs	80cm/31"	115cm/45"
		41222003	ClassicSling	polyester net	300kg/660lbs	XS	20–40kg/44–88lbs	30cm/12"
	41222004	ClassicSling	polyester net	300kg/660lbs	S	30–60kg/66–132lbs	38cm/15"	90cm/35"
	41222006	ClassicSling	polyester net	300kg/660lbs	M	55–85kg/121–187lbs	46cm/18"	95cm/37"
	41222007	ClassicSling	polyester net	300kg/660lbs	L	80–120kg/176–265lbs	63cm/25"	98cm/38"
	41222008	ClassicSling	polyester net	300kg/660lbs	XL	110–140kg/242–309lbs	80cm/31"	100cm/39"
		47001010	BariSling	polyester	570kg/1250lbs	XXXL	200–350kg/441–771lbs	–
47001011		BariSling	polyester	570kg/1250lbs	XXXXL	250–570kg/551–1257lbs	–	–
47091010		BariSling	disposable, non-woven	450kg/1000lbs	XXXL	200–450kg/441–1000lbs	–	–
	41221404	ClassicSling, with commode opening	polyester net	300kg/660lbs	S	30–60kg/66–132lbs	38cm/15"	90cm/35"
	41221406	ClassicSling, with commode opening	polyester net	300kg/660lbs	M	55–85kg/121–187lbs	46cm/18"	95cm/37"
	41221407	ClassicSling, with commode opening	polyester net	300kg/660lbs	L	80–120kg/132–198lbs	63cm/25"	98cm/38"





To and from toilet/High and shoulder-high back

Art. No.	Product name	Material	SWL	Size	Rec user weight	Waist width user	Sling overall height
43704002	HygieneHBSling	polyester	300kg/660lbs	XXS	10–30kg/22–66lbs	60cm/24"	92cm/36"
43704003	HygieneHBSling	polyester	300kg/660lbs	XS	25–45kg/55–99lbs	70cm/28"	92cm/36"
43704004	HygieneHBSling	polyester	300kg/660lbs	S	40–60kg/88–132lbs	80cm/31"	96cm/38"
43704006	HygieneHBSling	polyester	300kg/660lbs	M	55–80kg/121–176lbs	95cm/37"	102cm/40"
43704007	HygieneHBSling	polyester	300kg/660lbs	L	75–110kg/165–242lbs	110cm/43"	102cm/40"
43704008	HygieneHBSling	polyester	300kg/660lbs	XL	90–130kg/198–287lbs	130cm/51"	104cm/41"
43704009	HygieneHBSling	polyester	300kg/660lbs	XXL	120–170kg/265–376lbs	150cm/59"	104cm/41"
42521403	FlexibleSling, undivided, with commode opening	polyester net	300kg/660lbs	XS	20–40kg/44–88lbs	75cm/30"	80cm/31"
42521404	FlexibleSling, undivided, with commode opening	polyester net	300kg/660lbs	S	30–70kg/66–154lbs	85cm/33"	85cm/33"
42521406	FlexibleSling, undivided, with commode opening	polyester net	300kg/660lbs	M	60–90kg/132–198lbs	94cm/37"	93cm/37"
42521407	FlexibleSling, undivided, with commode opening	polyester net	300kg/660lbs	L	80–130kg/176–287lbs	105cm/41"	96cm/38"
42521408	FlexibleSling, undivided, with commode opening	polyester net	300kg/660lbs	XL	110–180kg/242–397lbs	115cm/45"	98cm/38"





To and from toilet/Low back

Art. No.	Product name	Material/Description	SWL	Size	Rec user weight	Sling seat measurement	Sling overall height
43504002	HygieneSling	polyester, with belt	300kg/660lbs	XXS	10–30kg/22–66lbs	25cm/10"	50cm/20"
43504003	HygieneSling	polyester, with belt	300kg/660lbs	XS	25–45kg/55–99lbs	35cm/14"	65cm/26"
43504004	HygieneSling	polyester, with belt	300kg/660lbs	S	40–60kg/88–132lbs	48cm/19"	75cm/30"
43504006	HygieneSling	polyester, with belt	300kg/660lbs	M	55–85kg/121–187lbs	55cm/22"	87cm/35"
43504007	HygieneSling	polyester, with belt	300kg/660lbs	L	75–100kg/165–220lbs	60cm/24"	93cm/37"
43504008	HygieneSling	polyester, with belt	300kg/660lbs	XL	90–130kg/198–287lbs	80cm/31"	95cm/37"
43504009	HygieneSling	polyester, with belt	300kg/660lbs	XXL	120–170kg/265–376lbs	90cm/35"	95cm/37"
43604002	HygieneSling	polyester, without belt	300kg/660lbs	XXS	10–30kg/22–66lbs	25cm/10"	50cm/20"
43604003	HygieneSling	polyester, without belt	300kg/660lbs	XS	25–45kg/55–99lbs	35cm/14"	65cm/26"
43604004	HygieneSling	polyester, without belt	300kg/660lbs	S	40–60kg/88–132lbs	48cm/19"	75cm/30"
43604006	HygieneSling	polyester, without belt	300kg/660lbs	M	55–85kg/121–187lbs	55cm/22"	87cm/35"
43604007	HygieneSling	polyester, without belt	300kg/660lbs	L	75–100kg/165–220lbs	60cm/24"	93cm/37"
43604008	HygieneSling	polyester, without belt	300kg/660lbs	XL	90–130kg/198–287lbs	80cm/31"	95cm/37"
43604009	HygieneSling	polyester, without belt	300kg/660lbs	XXL	120–170kg/265–376lbs	90cm/35"	95cm/37"
43904003	HygieneSling	polyester, extra back support	300kg/660lbs	XS	25–45kg/55–99lbs	35cm/14"	65cm/26"
43904004	HygieneSling	polyester, extra back support	300kg/660lbs	S	40–60kg/88–132lbs	48cm/19"	75cm/30"
43904006	HygieneSling	polyester, extra back support	300kg/660lbs	M	55–85kg/121–187lbs	55cm/22"	87cm/35"
43904007	HygieneSling	polyester, extra back support	300kg/660lbs	L	75–100kg/165–220lbs	60cm/24"	93cm/37"
43904008	HygieneSling	polyester, extra back support	300kg/660lbs	XL	90–130kg/198–287lbs	80cm/31"	95cm/37"
43804002	HygieneLBSling	polyester	300kg/660lbs	XXS	10–30kg/22–66lbs	60cm/24"	92cm/36"
43804003	HygieneLBSling	polyester	300kg/660lbs	XS	25–45kg/55–99lbs	70cm/28"	92cm/36"
43804004	HygieneLBSling	polyester	300kg/660lbs	S	40–60kg/88–132lbs	80cm/31"	96cm/38"
43804006	HygieneLBSling	polyester	300kg/660lbs	M	55–80kg/121–176lbs	95cm/37"	102cm/40"
43804007	HygieneLBSling	polyester	300kg/660lbs	L	75–110kg/165–242lbs	110cm/43"	102cm/40"
43804008	HygieneLBSling	polyester	300kg/660lbs	XL	90–130kg/198–287lbs	130cm/51"	104cm/41"
43804009	HygieneLBSling	polyester	300kg/660lbs	XXL	120–170kg/265–376lbs	150cm/59"	104cm/41"

Horizontal lifting/Turning

	Art. No.	Product name	Material/Description	SWL	Size	Rec user weight	Length	Width
	46502007	StretcherSling	polyester	300kg/660lbs	one size	50–200kg/110–441lbs	180cm/71"	86cm/34"
	46592007	StretcherSling	disposable, <i>non-woven</i>	300kg/660lbs	one size	50–200kg/110–441lbs	180cm/71"	86cm/34"
	47522007	SafeHandlingSheet	polyester net	300kg/660lbs	L	50–300kg/110–661lbs	240cm/94"	115cm/45"
	47542007	SafeHandlingSheet	polyester net/WendyLett4Way	300kg/660lbs	L	50–300kg/110–661lbs	240cm/94"	115cm/45"
	47502008	SafeHandlingSheet	polyester	400kg/900lbs	XL	50–400kg/110–881lbs	240cm/94"	115cm/45"
	47502009	SafeHandlingSheet	polyester, <i>reinforced</i>	570kg/1250lbs	XXL	50–570kg/110–1257lbs	240cm/94"	130cm/51"
	547810	PositioningSling	disposable, <i>non-woven, polyester</i>	454kg/1000lbs	one size	50–400kg/110–661lbs	198cm/78"	132cm/52"
	49592000	TurningSling, two loops	disposable, <i>non-woven</i>	250kg/550lbs	one size	50–250kg/110–551lbs	–	–
	49592005	TurningSling, one loop	disposable, <i>non-woven</i>	250kg/550lbs	one size	50–250kg/110–551lbs	–	–

Rising/Standing with MiniLift

	Art. No.	Product name	Material/Description	SWL	Size	Rec user weight	Lift vest width
	45500003	ThoraxSling	polyester	300kg/660lbs	XS	20–50kg/44–110lbs	85cm/33"
	45500004	ThoraxSling	polyester	300kg/660lbs	S	30–60kg/66–132lbs	94cm/37"
	45500006	ThoraxSling	polyester	300kg/660lbs	M	50–80kg/110–176lbs	104cm/41"
	45500007	ThoraxSling	polyester	300kg/660lbs	L	70–110kg/154–242lbs	124cm/49"
	45500008	ThoraxSling	polyester	300kg/660lbs	XL	100–150kg/220–330lbs	146cm/57"
	45500009	ThoraxSling	polyester	300kg/660lbs	XXL	140–200kg/309–441lbs	170cm/67"
	45550004	ThoraxSling	wipeable	205kg/450lbs	S	30–60kg/66–132lbs	94cm/37"
	45550006	ThoraxSling	wipeable	205kg/450lbs	M	50–80kg/110–176lbs	104cm/41"
	45550007	ThoraxSling	wipeable	205kg/450lbs	L	70–110kg/154–242lbs	124cm/49"
	45550008	ThoraxSling	wipeable	205kg/450lbs	XL	100–150kg/220–330lbs	146cm/57"
	45550009	ThoraxSling	wipeable	205kg/450lbs	XXL	140–200kg/309–441lbs	170cm/67"
	45600004	ThoraxSling	polyester, <i>with seat support</i>	300kg/660lbs	S	30–60kg/66–132lbs	94cm/37"
	45600006	ThoraxSling	polyester, <i>with seat support</i>	300kg/660lbs	M	50–80kg/110–176lbs	104cm/41"
	45600007	ThoraxSling	polyester, <i>with seat support</i>	300kg/660lbs	L	70–110kg/154–242lbs	124cm/49"
	45600008	ThoraxSling	polyester, <i>with seat support</i>	300kg/660lbs	XL	100–150kg/220–330lbs	146cm/57"
	45690004	ThoraxSling	disposable, <i>non-woven, with seat support</i>	200kg/440lbs	S	30–60kg/66–132lbs	94cm/37"
	45690006	ThoraxSling	disposable, <i>non-woven, with seat support</i>	200kg/440lbs	M	50–80kg/110–176lbs	104cm/41"
	45690007	ThoraxSling	disposable, <i>non-woven, with seat support</i>	200kg/440lbs	L	70–110kg/154–242lbs	124cm/49"
	45690008	ThoraxSling	disposable, <i>non-woven, with seat support</i>	200kg/440lbs	XL	100–150kg/220–330lbs	146cm/57"

Standing/Walking



Art No.	Product name	Material	SWL	Size	Rec user weight	Vest width	Lift vest overall height
44503003	StandingVest	polyester	300kg/660lbs	XS	15–40kg/33–88lbs	82cm/32"	68cm/27"
44503004	StandingVest	polyester	300kg/660lbs	S	20–50kg/44–110lbs	96cm/38"	73cm/29"
44503006	StandingVest	polyester	300kg/660lbs	M	40–60kg/88–132lbs	110cm/43"	73cm/29"
44503007	StandingVest	polyester	300kg/660lbs	L	50–80kg/110–176lbs	122cm/48"	73cm/29"
44503008	StandingVest	polyester	300kg/660lbs	XL	70–100kg/154–220lbs	130cm/51"	73cm/29"



44200002	WalkingVest	polyester	300kg/660lbs	XXS	15–40kg/33–88lbs	60cm/24"	67cm/26"
44200003	WalkingVest	polyester	300kg/660lbs	XS	25–50kg/55–110lbs	70cm/28"	68cm/27"
44200004	WalkingVest	polyester	300kg/660lbs	S	40–60kg/88–132lbs	82cm/32"	68cm/27"
44200006	WalkingVest	polyester	300kg/660lbs	M	55–80kg/121–176lbs	96cm/38"	73cm/29"
44200007	WalkingVest	polyester	300kg/660lbs	L	75–100kg/165–220lbs	110cm/43"	73cm/29"
44200008	WalkingVest	polyester	300kg/660lbs	XL	90–110kg/198–242lbs	122cm/48"	73cm/29"
44200009	WalkingVest	polyester	300kg/660lbs	XXL	100–130kg/220–287lbs	142cm/56"	75cm/30"



47103010	BariVest	polyester	570kg/1250lbs	XXXL	180–300kg/397–661lbs	190–280cm/75–110"	–
47103011	BariVest	polyester	570kg/1250lbs	XXXXL	250–570kg/551–1257lbs	230–310cm/91–122"	–

As a global supplier of equipment that is designed to facilitate better health comes, accelerate patient recovery and provide easier daily living for the elderly, disabled or those providing care and support to them, DHG is acutely aware that for equipment to be used to its full potential, those using it must have an understanding of everything a piece of equipment can do and how it should be used. The Ethos Education Programme underpins this theory by providing practical tips and techniques supported by CPD accredited theory-based learning for a wide range of healthcare practitioners.

Ethos webcasts

Sling Application, Part I:

Tips & Techniques for Sling Application In Bed

This Webcast focuses on the often-difficult challenge surrounding the application and assessment of slings. With an emphasis on best practice, this webcast explores the variety of solutions for sling application whilst your client is in bed, aiming to reduce the effort required to complete the task and ensure client comfort and dignity.

Sling Application, Part II:

Tips & Techniques for Sling Application In a Chair

This Webcast focuses on the often-difficult challenge surrounding the application and assessment of slings. With an emphasis on best practice, this webcast explores the variety of solutions for sling application whilst your client is seated, aiming to reduce the effort required to complete the task and ensure client comfort and dignity.

[Click here to view the Ethos Webcast Programme](#)



Withey Court, Western Industrial Estate
Caerphilly, United Kingdom CF83 1BF
T: +44 (0) 800 043 0881
F: +44 (0) 845 459 9832
info@directhealthcaregroup.com
DHG-HEALTHCARE.COM

SystemRoMedic™ is a complete solution for easy transfers.
The concept encompasses four product areas that meet
most needs in different transfer situations.

Transfer, between two locations. **Positioning**, within the same location.
Support, for mobility. **Lifting**, both manual and mechanical.

DHG
Moving Health Forward