

The logo for APEX, featuring the word "APEX" in a bold, sans-serif font. The letter "A" is stylized with a diagonal slash and a red dot above it. The background of the entire page is a low-angle shot of a modern building's facade with a grid of windows, set against a clear blue sky. The image is partially obscured by a dark teal vertical bar on the right side.

APEX


Pro-CARE 2

Pro-CARE 2

En idealisk madrass som förebygger och behandlar trycksår t.o.m sårkategori 2. Madrassen har special-designade rektangulära celler som ger ett bättre stöd för patienten och inbyggd hälavlastning i fotändan. Pumpen har ljud-och ljuslarm som indikerar om trycket blir för lågt i madrassen och/eller om det blir ström-avbrott. Viktindikatorn på panelen gör det enkelt för personalen att ställa in madrassens komfort.

■ **Alternerande läge**

Genom en periodisk 1-i-2-alternerande intervall med en tidsintervall på 15 minuter kan systemet undvika långvarigt tryck på lokal vävnad.

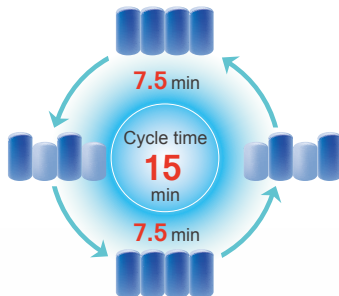

■ **Tryckjustering/komfortinställning**

Med viktindikator på panelen med ljud- och visuella larm, kan vårdgivare problemfritt ställa in en optimerad tryckinställning.

■ **Integrerad hälavlastning**

De fem nedersta cellerna i madrassen kan individuellt kopplas ur, beroende på patientens längd, för att helt eliminera trycket på hälsarna.



■ **Statiskt läge**

Kan användas i omvårdnadssituationer såsom lägesändringar eller vid omvårdnad. Återgår automatiskt till växlande läge efter 20 minuter.

■ **Rektangulära celler**

Cellerna anpassar sig efter kroppens konturer och rörelse för att minimera skjuvkrifter. Optimerar tryckfördelningen.


Specifikationer

	Modell	Pro-CARE 2
Pump	Storlek	29,1 x 20 x 11,7 cm
	Vikt	2,2 Kg
	Höljets material	Brandhämmande ABS
	Matningsspänning	AC 220-240V, 50Hz; AC 120V/ 60Hz
	Operationsintervall	15 minuter
Madrass	Särkategori	2
	Storlek	200x80/85/90x12,7cm
	Celler	18 x 12,5 cm höga celler
	Vikt	7,3 kg
	Överdrag	2-vägs stretch PU
	Cellmaterial	TPU
	Brukarvikt	20-180 kg


Pro-CARE 2

Subjekt: 176 cm/82 kg

Tryckkarta 1

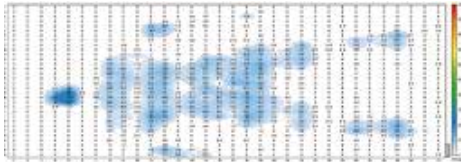

Tryckkarta 2

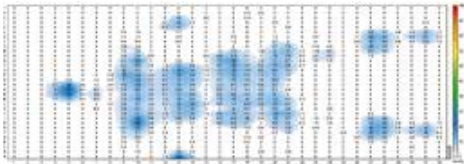

Pyo-CARE series Therapeutic Support Surface



APEX MEDICAL CORP.

www.apexmedicalcorp.com

HEADQUARTERS

APEX MEDICAL CORP.

No. 9, Min Sheng St., Tu-Cheng, New Taipei City, 23679, Taiwan

TEL:+886 2 2268 3100

FAX:+886 2 2268 6525

marketing@apexmedical.com.tw

EUROPE

APEX MEDICAL S.L.

Elcano 9, 6ª planta 48008 Bilbao, Vizcaya, Spain

TEL:+34 94 4706408

FAX:+34 94 4706409

info@apexmedical.es

UK

APEX MEDICAL LIMITED

Atlantic House, Imperial Way Reading

Berkshire RG2 0TD, United Kingdom

TEL:+44(0)118 9908 049

FAX:+44(0)118 9036 100

sales@apexmedical.org.uk

MIDDLE EAST • AFRICA • INDIA

29, Misr Liltameer Building, Stree#10, AREA#8, Sheraton Building,

Heliopolis Cairo, Egypt.

TEL:+202 2 2671303

FAX:+202 2 2678649

sunmed@link.net

USA

APEX MEDICAL USA CORP.

615 N. Berry St. Suite D, Brea, CA 92821, USA

TEL: +1 714 671 3818

FAX: +1 714 671 9040

sales.usa@apexmedicalcorp.com

CHINA

APEX MEDICAL (KUNSHAN) CORP.

No. 1368, ZiZhu Road, Kunshan,

Jiang Su Province, China

TEL:+86 512 8617 1660 Ext1150

FAX:+86 512 8617 8663

customer@apexmedical.com.cn

Återförsäljare

Vill du veta mer?

Kontakta oss på: tel 0371-318 00

Gate Rehab Development AB

Industrigatan 2 | S-333 75 Reftele

info@gaterd.se | www.gaterd.se

GATE

WE DO BETTER.