MiniLift160EM/EE





A unique solution for active sit-to-stand - four variants available

MiniLift160 is a mobile sit-to-stand lift which is used, in combination with a specially adapted lifting vest, to assist users with reduced function when rising from a sitting to a standing position. MiniLift160 features a construction offering a totally unique and extremely comfortable sit-to-stand experience. The lift is also compact and convenient and very easy to manoeuvre; features which, in combination with the possibility to adjust the settings for an optimal fit for every user, makes it very useful in most locations and situations. MiniLift160 offers a unique combination of functionality, performance and quality in relation to price. A more economical solution is hard to find!





MiniLift160EM/EE

Extremely comfortable, active sit-to-stand with natural movement pattern

MiniLift160 has been developed to, as gently and comfortably as possible, assist the user when rising from a sitting to a standing position. The handle bar and the side supports offer multiple hand placement options for the user. The stable anti-slip foot plate is very low and the lower leg supports, as well as the lifting arm, can be adjusted to fit every user optimally. When MiniLift160 is used in combination with the appropriate lifting accessories, the user is supported under the feet, for the front of the lower legs and behind the back, which provides for a safe and secure, yet active sit-to-stand procedure. The construction safely moves the user forwards and upwards in a natural movement pattern and, at the same time, the leg muscles and balance are exercised. MiniLift160 is available in four variants with electrical or manual base widening and with standard height or low legs for use with lower beds. All variants are manufactured in steel for stability and strength and are suitable for users weighing up to 160 kg/350 lbs.





New control box with diagnostics, e.g. lift counter and overload

indicator

Electrical emergency stop button readily available on the control hox


PRODUCT SPECIFICATION

| Product name Description | Art. No. |
|--|--------------------|
| MiniLift160EMMobile sit-to-stand lift with manual base widening | 60300010 |
| MiniLift160EMMobile sit-to-stand lift with manual base widening, low leg | js <u>60300011</u> |
| MiniLift160EE Mobile sit-to-stand lift with electrical base widening | 60300012 |
| MiniLift160EE Mobile sit-to-stand lift with electrical base widening, low le | gs 60300013 |

TECHNICAL INFORMATION

| Safe working load (SWL) | 160 kg / 350 lbs |
|--|---------------------------|
| User length, min-max | 1400-2000 mm / 55.1-78.7 |
| Material | Powder-coated stee |
| Height | 980 mm / 38.6 |
| Length | 919 mm / 36.2 |
| Length, folded | 1098 mm / 43.2 |
| Width of base, outer measurement | 669-959 mm / 26.3-37.8 |
| Width of base, inner measurement | 531-773 mm / 20.9-30.4 |
| Weight | 40.5 kg / 89.3 lbs |
| Base widening, EM | Manua |
| Base widening, EE | Electrica |
| Height of footplate, for user | 59 mm / 2.3 |
| Castors, diameter, 60300010 and 60300012 | 100 mm / 3.9 [°] |
| Rear castors, diameter, 60300011 and 60300013 | 100 mm / 3.9 [°] |
| Front castors, diameter, 60300011 and 60300013 | 70 mm / 2.8 |
| Lifting speed, without load | 35 mm / 1.4" per sec |
| Emergency lowering | Manual and electrica |

ACCESSORIES

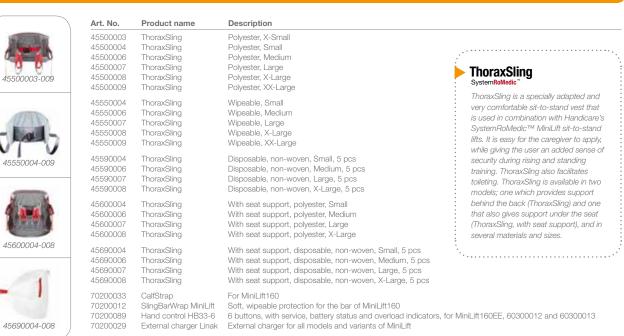

Other models and variants of MiniLift mobile sit-to-stand lifts

MiniLift is available in two more models with different lifting capacity. Both feature electrical base-width adjustment and standard height legs.

| Art. No. | Product name |
|-----------|--------------|
| 401100334 | MiniLift200 |
| 400641434 | MiniLift125 |

Description Mobile sit-to-stand lift with electrical base widening, SWL 200 kg/440 lbs Mobile sit-to-stand lift with electrical base widening, SWL 125 kg/275 lbs

Please visit www.handicare.co.uk for more information.



handicare 82 First avenue Pensnett Estate. Kingswinford West Midlands DY6 7E.I T: 01384 408700 www.handicare.co.uk enquiries@handicare.co.uk

SystemRoMedic™ is a complete solution for easy transfers. The concept encompasses four product areas that meet most needs in different transfer situations. Transfer, between two locations. Positioning, within the same location. Support, for mobility. Lifting, both manual and mechanical.