

RisePorto300™

SystemRoMedic™



SWL
140 kg
300 lbs

Safe and flexible lifting solutions for all settings and situations

RisePorto300 is a portable ceiling lift unit featuring unusually small measurements and extremely low weight. At the same time, it is both strong and safe, providing secure, comfortable and individually adapted lifting solutions for all situations and for users weighing up to 140 kg/300 lbs. RisePorto300 can be used on all rail systems, for all settings and in all types of rooms. The lift unit can easily be moved between different rail systems and rooms. This supports a high degree of utilization which, in turn, avoids manual lifting and work-related injuries and accidents. RisePorto300 is a modern and stylish ceiling lift unit offering safe and flexible lifting solutions characterized by performance, simplicity and comfort. It also contributes to the creation of professional and highly functional care environments with a harmonic and human atmosphere.



PRODUCT SPECIFICATION

| Item No. | Product | Description |
|----------|--------------|--------------------------------------------------------------|
| 50200016 | RisePorto300 | Complete with hand control and charger, SWL 140 kg / 300 lbs |

Trolley ordered separately.


ACCESSORIES

| Item No. | Product | Description |
|----------|--------------------------|----------------------|
| 50200014 | Trolley 360 RisePorto | 360° turning trolley |
| 50200004 | Room transfer strap | 300 mm / 11.8" |
| 70200021 | Room transfer strap | 700 mm / 27.6" |
| 70200022 | Room transfer strap | 800 mm / 31.5" |
| 70200023 | Room transfer strap | 900 mm / 35.4" |
| 70200024 | Room transfer strap | 1000 mm / 39.4" |
| 50200017 | Handle for mounting | 750 mm / 29.5" |
| 70200053 | Lift strap extension | 370 mm / 14.6" |
| 70200054 | Lift strap extension | 570 mm / 22.4" |
| 70200055 | Lift strap extension | 670 mm / 26.4" |
| 70200056 | Lift strap extension | 770 mm / 30.3" |
| 70200057 | Lift strap extension | 870 mm / 34.2" |
| 70200058 | Lift strap extension | 970 mm / 38.2" |
| 70200048 | Software for diagnostics | Software and cable |

TECHNICAL INFORMATION

| | |
|----------------------------------------|---------------------------------------|
| Safe working load, RisePorto300: | 140 kg / 300 lbs |
| Material, chassis: | Aluminium |
| Material, cover: | Non flammable ABS plastic |
| Material, sling bar: | Aluminium |
| Measurements, length x width x height: | 288 x 152 x 153 mm / 11.3 x 5.98 x 6" |
| Build in height: | 270 mm / 10.6" |
| Weight: | 5.8 kg / 12.8 lbs |
| Transfer in rail: | Manual |
| Sling bar: | Integrated, foldable |
| Sling bar, width: | 458 mm / 18" |
| Hand control: | Electrical, 4 buttons |
| Lifting range: | 1700 mm / 66.9" |
| Lifting speed, without load: | 30 mm / 1.2" per sec |
| Lifting speed, with load: | 23 mm / 0.9" per sec |
| Sound level, with load: | 57 dB(A) |
| Emergency lowering: | Electrical and manual |
| Emergency stop: | Manual |
| Battery: | NiMH 26.4V 2.7 Ah |
| Protection class, lift unit: | IPX3 |
| Protection class, hand control: | IPX4 |

CE RisePorto300 meets the requirements of the Medical Devices Directive for Class 1 products 93/42/EEC.

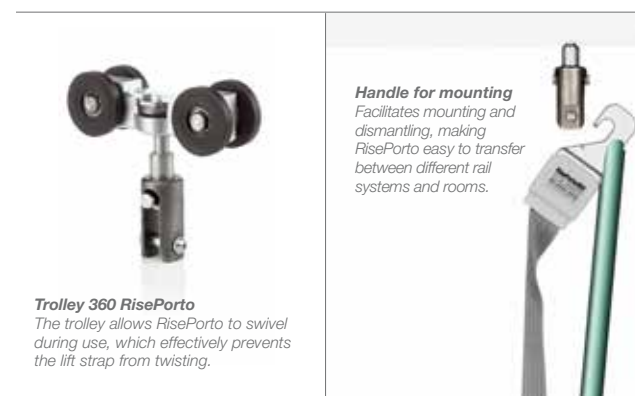

Trolley 360 RisePorto
The trolley allows RisePorto to swivel during use, which effectively prevents the lift strap from twisting.

Handle for mounting
Facilitates mounting and dismantling, making RisePorto easy to transfer between different rail systems and rooms.

Please visit www.handicare.co.uk for more information.

SystemRoMedic™ is a complete solution for easy transfers. The concept encompasses four product areas that meet most needs in different transfer situations.

Transfer, between two locations. **Positioning**, within the same location.

Support, for mobility. **Lifting**, both manual and mechanical.



handicare

82 First avenue
Pensnett Estate,
Kingswinford
West Midlands DY6 7FJ
T: 01384 408700
www.handicare.co.uk
enquiries@handicare.co.uk