

RisePorto450™

SystemRoMedic™

SWL
205 kg
450 lbs



Safe and flexible lifting solutions for all settings and situations

RisePorto450 is a portable ceiling lift unit featuring unusually small measurements and extremely low weight. At the same time, it is both strong and safe, providing secure, comfortable and individually adapted lifting solutions for all situations and for users weighing up to 205 kg/450 lbs. RisePorto450 can be used on all rail systems, for all settings and in all types of rooms. The lift unit can easily be moved between different rail systems and rooms. This supports a high degree of utilization which, in turn, avoids manual lifting and work-related injuries and accidents. RisePorto450 is a modern and stylish ceiling lift unit offering safe and flexible lifting solutions characterized by performance, simplicity and comfort. It also contributes to the creation of professional and highly functional care environments with a harmonic and human atmosphere.



RisePorto450™

PRODUCT SPECIFICATION

Item No.	Product	Description
50200015	RisePorto450	Complete with hand control and charger, SWL 205 kg / 450 lbs

Trolley ordered separately.


ACCESSORIES

Item No.	Product	Description
50200014	Trolley 360 RisePorto	360° turning trolley
50200004	Room transfer strap	300 mm / 11.8"
70200021	Room transfer strap	700 mm / 27.6"
70200022	Room transfer strap	800 mm / 31.5"
70200023	Room transfer strap	900 mm / 35.4"
70200024	Room transfer strap	1000 mm / 39.4"
50200017	Handle for mounting	750 mm / 29.5"
70200053	Lift strap extension	370 mm / 14.6"
70200054	Lift strap extension	570 mm / 22.4"
70200055	Lift strap extension	670 mm / 26.4"
70200056	Lift strap extension	770 mm / 30.3"
70200057	Lift strap extension	870 mm / 34.2"
70200058	Lift strap extension	970 mm / 38.2"
70200048	Software for diagnostics	Software and cable

TECHNICAL INFORMATION

Safe working load, RisePorto450:	205 kg / 450 lbs
Material, chassis:	Aluminium
Material, cover:	Non flammable ABS plastic
Material, sling bar:	Aluminium
Measurements, length x width x height:	288 x 152 x 153 mm / 11.3 x 5.98 x 6"
Build in height:	270 mm / 10.6"
Weight:	5.8 kg / 12.8 lbs
Transfer in rail:	Manual
Sling bar:	Integrated, foldable
Sling bar, width:	458 mm / 18"
Hand control:	Electrical, 4 buttons
Lifting range:	1700 mm / 66.9"
Lifting speed, without load:	30 mm / 1.2" per sec
Lifting speed, with load:	18 mm / 0.7" per sec
Sound level, with load:	57 dB(A)
Emergency lowering:	Electrical and manual
Emergency stop:	Manual
Battery:	NiMH 26.4V 2.7 Ah
Protection class, lift unit:	IPX3
Protection class, hand control:	IPX4

CE RisePorto450 meets the requirements of the Medical Devices Directive for Class 1 products 93/42/EEC.

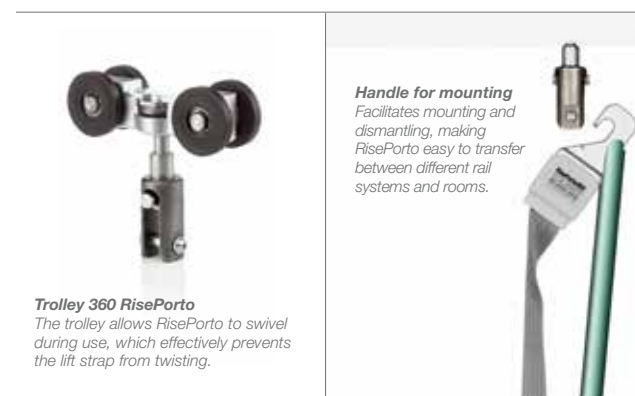

Trolley 360 RisePorto
The trolley allows RisePorto to swivel during use, which effectively prevents the lift strap from twisting.

Handle for mounting
Facilitates mounting and dismantling, making RisePorto easy to transfer between different rail systems and rooms.

Please visit www.handicare.co.uk for more information.



handicare

82 First avenue
Pensnett Estate,
Kingswinford
West Midlands DY6 7FJ
T: 01384 408700
www.handicare.co.uk
enquiries@handicare.co.uk

SystemRoMedic™ is a complete solution for easy transfers. The concept encompasses four product areas that meet most needs in different transfer situations.

Transfer, between two locations. **Positioning**, within the same location.

Support, for mobility. **Lifting**, both manual and mechanical.